

Product Datasheet

Lightning-Link (R) Rapid DyLight 594 Antibody Labeling Kit

324-0015

Unit Size: 1 mg

Store at -20C.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/324-0015

Updated 5/23/2024 v.20.1

Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/324-0015



324-0015**Lightning-Link (R) Rapid DyLight 594 Antibody Labeling Kit**


Product Information	
Unit Size	1 mg
Concentration	Concentration is not relevant for this product. Please see the protocols for proper use of this product.
Storage	Store at -20C.
Conjugate	DyLight 594
Product Description	
Description	<p>Lightning-Link Rapid is an innovative technology that enables direct labeling of proteins, peptides or other biomolecules for use in R&D applications, drug discovery and the development of diagnostic kits (See protocol for further information).</p> <p>The easy-to-use, one step procedure allows researchers to covalently label biomolecules with only 30 seconds hands-on time; furthermore conjugates are ready to use in less than twenty minutes.</p> <p>The researcher simply pipettes the biomolecule into a vial of lyophilized mixture containing the label of interest and incubates (for more details please watch the video below).</p> <p>Features Quick and easy to use Benefits Save time, no special knowledge required No separation steps 100% recovery - no antibody/protein loss Can be used in a wide range of applications Flexible Freeze dried Ships at ambient temperature, long shelf-life Fully scalable (10 ug to 1 g or more) Easy transfer from R&D to manufacturing Stringently QC tested Consistent high quality, excellent batch-to-batch reproducibility Large number of labels available Experimental flexibility Reliable: nearly 300 references Successfully used in many fields of research</p> <p>DyLight 594 generates red fluorescence for a wide array of fluorescence labeling-based applications. It has a strong absorption at 594nm, high fluorescence at 629nm (extinction coefficient $8.0 \times 10^4 \text{ cm}^{-1}\text{M}^{-1}$) and high quantum yield.</p> <p>Learn more about Lightning-Link™ Conjugation Kits by reading FAQs</p> <p>For more information please check out these useful links! Antibody Labeling Guide Antibody Conjugation Illustrated Assay</p>
Kit Components	1 or 3 glass vial(s) of Lightning-Link Rapid mix, 1 vial of LL-Rapid Modifier reagent, 1 vial of LL-Rapid Quencher reagent
Notes	<p>This product is manufactured by Abcam and distributed by Novus Biologicals.</p> <p>This product is for research use only and is not approved for use in humans or in clinical diagnosis. This product is guaranteed for 1 year from date of receipt and this statement overrides any mentioned guarantee period on the limitations section of this products datasheet. Please contact technical@novusbio.com with questions.</p>

Product Application Details


Applications	Flow Cytometry
Recommended Dilutions	Flow Cytometry
Application Notes	By circumventing the desalting or dialysis steps that commonly interrupt traditional antibody conjugation procedures, LightningLink technology can be used to label both small (e.g. 10 ug) and large quantities of primary antibodies with ease. Batch-to-batch variation upon scale up is minimal as the process is so simple, and recoveries are always 100%.

Images

Flow Cytometry: Lightning-Link Rapid DyLight 594 Antibody Labeling Kit [324-0015] - Mouse anti-human CD3 was conjugated with DyLight® 594 using Lightning-Link® Rapid kit. The conjugated antibody was then used to stain human peripheral blood lymphocytes, followed by analysis with flow cytometry. (Blue line - negative control; red line - positive staining).

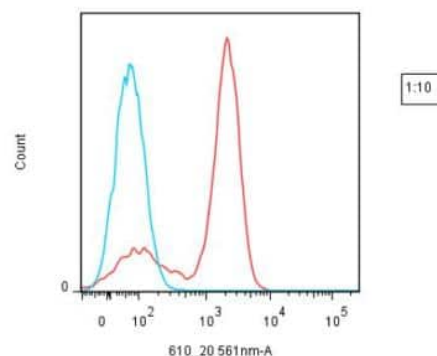

Lightning-Link Rapid DyLight 594 Antibody Labeling Kit [324-0015]

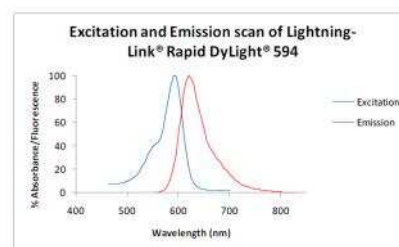

Scan performed in TBS pH 8.0

Absorbance Max (nm)	Emission Max (nm)	Extinction Coefficient (cm ² M ⁻¹)	Fluorescent Colour	Stokes Shift
593	618	80000	Red	25



Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112
USA
Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave
Toronto, ON M8Z 4E6
Canada
Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB, United Kingdom
Phone: (44) (0) 1235 529449
Free Phone: 0800 37 34 15
Fax: (44) (0) 1235 533420
info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com
Technical Support: nb-technical@bio-techne.com
Orders: nb-customerservice@bio-techne.com
General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Kits are guaranteed for 6 months from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/324-0015

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

