

Product Datasheet

Lightning-Link (R) PerCP Antibody Labeling Kit 718-0010

Unit Size: 3 x 100ug Reaction

Store at -20C.

www.novusbio.com



technical@novusbio.com

Publications: 9

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/718-0010

Updated 7/23/2023 v.20.1

**Earn rewards for product
reviews and publications.**

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/718-0010



718-0010**Lightning-Link (R) PerCP Antibody Labeling Kit**

Product Information	
Unit Size	3 x 100ug Reaction
Concentration	Concentration is not relevant for this product. Please see the protocols for proper use of this product.
Storage	Store at -20C.
Conjugate	PerCP

Product Description	
Description	<p>Lightning-Link antibody labeling kits enable the direct labeling of antibodies, proteins, peptides or other biomolecules for use in R&D applications, drug discovery and the development of diagnostic kits (See protocol for further information).</p> <p>Our PerCP antibody labeling kit enables the direct conjugation of PerCP to any biomolecule with an available amine group. The researcher simply pipettes the antibody or other biomolecule into the vial of Lightning-Link PerCP and incubates for 3 hours.</p> <p>FeaturesQuick and easy to use BenefitsSave time, no special knowledge required No separation steps 100% recovery - no antibody/protein loss Can be used in a wide range of applications Flexible Freeze dried Ships at ambient temperature, long shelf-life Fully scalable (10 ug to 1 g or more) Easy transfer from R&D to manufacturing Stringently QC tested Consistent high quality, excellent batch-to-batch reproducibility Large number of labels available Experimental flexibility Reliable: nearly 300 references Successfully used in many fields of research</p> <p>PerCP is a small 35kDa phycobiliprotein isolated from red algae. It has a strong absorption peak at 484nm and a large stokes shift (194nm) with its maximum emission peak at 678nm. PerCP is commonly used for fluorescent immunolabeling, particularly in applications involving fluorescent-activated cell sorting.</p> <p>Learn more about Lightning-Link™ Conjugation Kits by reading FAQs</p> <p>For more information please check out these useful links! Antibody Labeling Guide Antibody Conjugation Illustrated Assay</p>
Kit Components	1 or 3 glass vial(s) of Lightning-Link mix, 1 vial of LL-Modifier reagent, 1 vial of LL-Quencher reagent
Notes	<p>This product is manufactured by Abcam and distributed by Novus Biologicals.</p> <p>This product is for research use only and is not approved for use in humans or in clinical diagnosis. This product is guaranteed for 1 year from date of receipt and this statement overrides any mentioned guarantee period on the limitations section of this products datasheet. Please contact technical@novusbio.com with questions.</p>

Product Application Details	
Applications	Flow Cytometry
Recommended Dilutions	Flow Cytometry

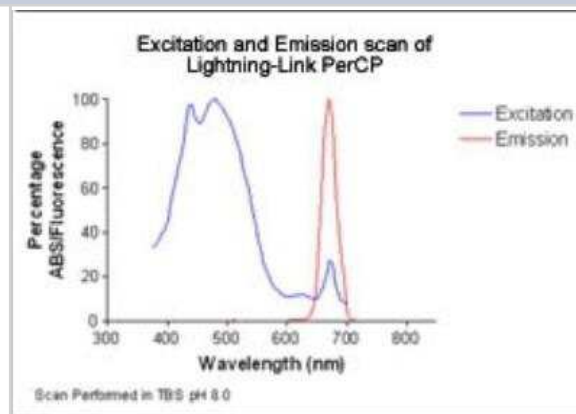


Application Notes

By circumventing the desalting or dialysis steps that commonly interrupt traditional antibody conjugation procedures, LightningLink technology can be used to label both small (e.g. 10 ug) and large quantities of primary antibodies with ease. Batch-to-batch variation upon scale up is minimal as the process is so simple, and recoveries are always 100%. This kit is supplied with 3 vials, each suitable for labeling up to 100 ug of antibody.

Images

Flow Cytometry: Lightning-Link PerCP Antibody Labeling Kit [718-0010]

**Publications**

Keeley EC, Schutt RC, Marinescu MA et al. Circulating fibrocytes as predictors of adverse events in unstable angina Transl Res. 2016-03-08 [PMID: 27012475] (FLOW)

Starkie DO, Compson JE, Rapecki S, Lightwood DJ. Generation of Recombinant Monoclonal Antibodies from Immunised Mice and Rabbits via Flow Cytometry and Sorting of Antigen-Specific IgG+ Memory B Cells PLoS One 2016-03-29 [PMID: 27022949] (FLOW)

Rodgers JM, Robinson AP, Rosler ES et al. IL-17A activates ERK1/2 and enhances differentiation of oligodendrocyte progenitor cells. Glia 2014-01-01 [PMID: 25557204]

Ipp H, Nkambule BB, Reid TD et al. CD4+ T cells in HIV infection show increased levels of expression of a receptor for vasoactive intestinal peptide, VPAC2. Immunol Res 2014-01-01 [PMID: 24469917] (FLOW)

Knudsen T, Kjalke M, Tranholm M et al. Development of a flow cytometric assay for detection of coated platelets in dogs and evaluation of binding of coated platelets to recombinant human coagulation factor VIIa American Journal of Veterinary Research 2011-01-01 [PMID: 21801056] (FLOW)

Al-Dujaili EA, Mullins LJ, Bailey MA, Kenyon CJ. Development of a Highly Sensitive ELISA for Aldosterone in Mouse Urine: Validation in Physiological Pathophysiological States of Aldosterone Excess Depletion - . Steroids. 2008-01-01 [PMID: 19162057]

Bao S, Wu Q, Li Z et al. Targeting Cancer Stem Cells through L1CAM Suppresses Glioma Growth. Cancer Res 2008-01-01 [PMID: 18676824]

Marr AK, Jenssen H, Moniri MR et al. Bovine lactoferrin and lactoferricin interfere with intracellular trafficking of Herpes simplex virus-1. Biochemie 2009-01-01 [PMID: 18573311]

Velappan N, Clements J, Kiss C et al. Fluorescence linked immunosorbant assays using microtiter plates. J Immunol Methods. 2013-01-01 [PMID: 18514691]



Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112

USA

Phone: 303.730.1950

Toll Free: 1.888.506.6887

Fax: 303.730.1966

nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave

Toronto, ON M8Z 4E6

Canada

Phone: 905.827.6400

Toll Free: 855.668.8722

Fax: 905.827.6402

canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane

Abingdon Science Park

Abingdon, OX14 3NB, United Kingdom

Phone: (44) (0) 1235 529449

Free Phone: 0800 37 34 15

Fax: (44) (0) 1235 533420

info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com

Technical Support: nb-technical@bio-techne.com

Orders: nb-customerservice@bio-techne.com

General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Kits are guaranteed for 6 months from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/718-0010

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

