

Product Datasheet

Lightning-Link (R) PerCP Antibody Labeling Kit 718-0030

Unit Size: 3 x 10ug Reaction

Store at -20C.

www.novusbio.com



technical@novusbio.com

Publications: 5

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/718-0030

Updated 10/23/2024 v.20.1

Earn rewards for product
reviews and publications.

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/718-0030



718-0030**Lightning-Link (R) PerCP Antibody Labeling Kit**

Product Information	
Unit Size	3 x 10ug Reaction
Concentration	Concentration is not relevant for this product. Please see the protocols for proper use of this product.
Storage	Store at -20C.
Conjugate	PerCP

Product Description	
Description	<p>Lightning-Link antibody labeling kits enable the direct labeling of antibodies, proteins, peptides or other biomolecules for use in R&D applications, drug discovery and the development of diagnostic kits (See protocol for further information).</p> <p>Our PerCP antibody labeling kit enables the direct conjugation of PerCP to any biomolecule with an available amine group. The researcher simply pipettes the antibody or other biomolecule into the vial of Lightning-Link PerCP and incubates for 3 hours.</p> <p>FeaturesQuick and easy to use BenefitsSave time, no special knowledge required No separation steps 100% recovery - no antibody/protein loss Can be used in a wide range of applications Flexible Freeze dried Ships at ambient temperature, long shelf-life Fully scalable (10 ug to 1 g or more) Easy transfer from R&D to manufacturing Stringently QC tested Consistent high quality, excellent batch-to-batch reproducibility Large number of labels available Experimental flexibility Reliable: nearly 300 references Successfully used in many fields of research</p> <p>PerCP is a small 35kDa phycobiliprotein isolated from red algae. It has a strong absorption peak at 484nm and a large stokes shift (194nm) with its maximum emission peak at 678nm. PerCP is commonly used for fluorescent immunolabeling, particularly in applications involving fluorescent-activated cell sorting.</p> <p>Learn more about Lightning-Link™ Conjugation Kits by reading FAQs</p> <p>For more information please check out these useful links! Antibody Labeling Guide Antibody Conjugation Illustrated Assay</p>
Kit Components	1 or 3 glass vial(s) of Lightning-Link mix, 1 vial of LL-Modifier reagent, 1 vial of LL-Quencher reagent
Notes	<p>This product is manufactured by Abcam and distributed by Novus Biologicals.</p> <p>This product is for research use only and is not approved for use in humans or in clinical diagnosis. This product is guaranteed for 1 year from date of receipt and this statement overrides any mentioned guarantee period on the limitations section of this products datasheet. Please contact technical@novusbio.com with questions.</p>

Product Application Details	
Applications	Flow Cytometry
Recommended Dilutions	Flow Cytometry

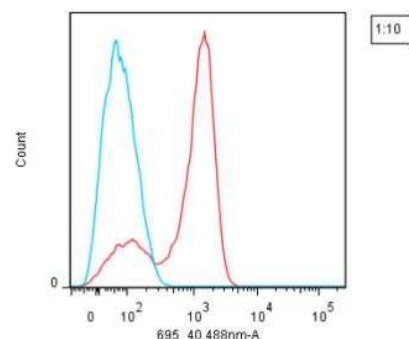


Application Notes

By circumventing the desalting or dialysis steps that commonly interrupt traditional antibody conjugation procedures, LightningLink technology can be used to label both small (e.g. 10 ug) and large quantities of primary antibodies with ease. Batch-to-batch variation upon scale up is minimal as the process is so simple, and recoveries are always 100%. This kit is supplied with 3 vials, each suitable for labeling up to 10 ug of antibody.

Images

Flow Cytometry: Lightning-Link PerCP Antibody Labeling Kit [718-0030] - Mouse anti-human CD3 was conjugated with PerCP using a Lightning-Link kit. The conjugated antibody was then used to stain human peripheral blood lymphocytes, followed by analysis with flow cytometry. (Blue line - negative control; red line - positive staining).

**Publications**

Starkie DO, Compson JE, Rapecki S, Lightwood DJ. Generation of Recombinant Monoclonal Antibodies from Immunised Mice and Rabbits via Flow Cytometry and Sorting of Antigen-Specific IgG+ Memory B Cells PLoS One 2016-03-29 [PMID: 27022949] (FLOW)

Keeley EC, Schutt RC, Marinescu MA et al. Circulating fibrocytes as predictors of adverse events in unstable angina Transl Res. 2016-03-08 [PMID: 27012475] (FLOW)

Rodgers JM, Robinson AP, Rosler ES et al. IL-17A activates ERK1/2 and enhances differentiation of oligodendrocyte progenitor cells. Glia 2014-01-01 [PMID: 25557204]

Ipp H, Nkambule BB, Reid TD et al. CD4+ T cells in HIV infection show increased levels of expression of a receptor for vasoactive intestinal peptide, VPAC2. Immunol Res 2014-01-01 [PMID: 24469917] (FLOW)

Knudsen T, Kjalke M, Tranholm M et al. Development of a flow cytometric assay for detection of coated platelets in dogs and evaluation of binding of coated platelets to recombinant human coagulation factor VIIa American Journal of Veterinary Research 2011-01-01 [PMID: 21801056] (FLOW)



Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112
USA
Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave
Toronto, ON M8Z 4E6
Canada
Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB, United Kingdom
Phone: (44) (0) 1235 529449
Free Phone: 0800 37 34 15
Fax: (44) (0) 1235 533420
info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com
Technical Support: nb-technical@bio-techne.com
Orders: nb-customerservice@bio-techne.com
General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Kits are guaranteed for 6 months from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/718-0030

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

