

# Product Datasheet

## Lightning-Link (R) APC-Cy7 Antibody Labeling Kit 765-0015

Unit Size: 1 mg

Store at -20C.

[www.novusbio.com](http://www.novusbio.com)



[technical@novusbio.com](mailto:technical@novusbio.com)

### Publications: 6

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:  
[www.novusbio.com/765-0015](http://www.novusbio.com/765-0015)

Updated 7/23/2023 v.20.1

Earn rewards for product  
reviews and publications.

Submit a publication at [www.novusbio.com/publications](http://www.novusbio.com/publications)

Submit a review at [www.novusbio.com/reviews/destination/765-0015](http://www.novusbio.com/reviews/destination/765-0015)



**765-0015****Lightning-Link (R) APC-Cy7 Antibody Labeling Kit**

<b>Product Information</b>	
<b>Unit Size</b>	1 mg
<b>Concentration</b>	Concentration is not relevant for this product. Please see the protocols for proper use of this product.
<b>Storage</b>	Store at -20C.
<b>Conjugate</b>	Allophycocyanin/Cy7
<b>Product Description</b>	
<b>Description</b>	<p>Lightning-Link antibody labeling kits enable the direct labeling of antibodies, proteins, peptides or other biomolecules for use in R&amp;D applications, drug discovery and the development of diagnostic kits (See protocol for further information).</p> <p>The easy-to-use, one step procedure allows researchers to covalently label biomolecules with only 30 seconds hands on time. The researcher simply pipettes the biomolecule into a vial of lyophilized mixture containing the label of interest and incubates.</p> <p><b>Features</b>Quick and easy to use  <b>Benefits</b>Save time, no special knowledge required            No separation steps            100% recovery - no antibody/protein loss            Can be used in a wide range of applications            Flexible            Freeze dried            Ships at ambient temperature, long shelf-life            Fully scalable (10 ug to 1 g or more)            Easy transfer from R&amp;D to manufacturing            Stringently QC tested            Consistent high quality, excellent batch-to-batch reproducibility            Large number of labels available            Experimental flexibility            Reliable: nearly 300 references            Successfully used in many fields of research</p> <p>APC/Cy7 is a tandem conjugate. The APC is excited at 652nm and functions as an energy donor for the Cy7. Energy is transferred from the APC to the Cy7 via <b>FRET</b> (Fluorescence Resonance Energy Transfer). The Cy7 emits the energy received from the APC in the form of long wavelength light at 790nm</p> <p>Learn more about Lightning-Link™ Conjugation Kits by reading <a href="#">FAQs</a></p> <p>For more information please check out these useful links!  <a href="#">Antibody Labeling Guide</a>  <a href="#">Antibody Conjugation Illustrated Assay</a></p>
<b>Kit Components</b>	1 or 3 glass vial(s) of Lightning-Link mix, 1 vial of LL-Modifier reagent, 1 vial of LL-Quencher reagent
<b>Notes</b>	<p>This product is manufactured by Abcam and distributed by Novus Biologicals.</p> <p>This product is for research use only and is not approved for use in humans or in clinical diagnosis. This product is guaranteed for 1 year from date of receipt and this statement overrides any mentioned guarantee period on the limitations section of this products datasheet. Please contact <a href="mailto:technical@novusbio.com">technical@novusbio.com</a> with questions.</p>
<b>Product Application Details</b>	
<b>Applications</b>	Flow Cytometry
<b>Recommended Dilutions</b>	Flow Cytometry

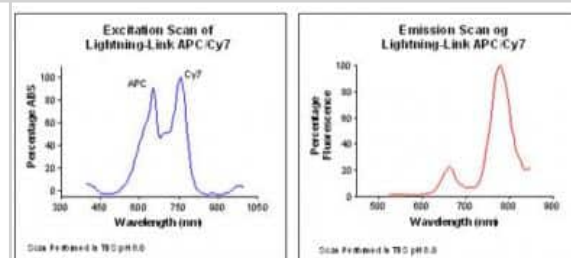


**Application Notes**

By circumventing the desalting or dialysis steps that commonly interrupt traditional antibody conjugation procedures, LightningLink technology can be used to label both small (e.g. 10 ug) and large quantities of primary antibodies with ease. Batch-to-batch variation upon scale up is minimal as the process is so simple, and recoveries are always 100%. This kit can be used to label up to 1.5 mg of antibody, and is supplied in one vial.

**Images**

Flow Cytometry: Lightning-Link APC-Cy7 Antibody Labeling Kit [765-0015]

**Publications**

Okagawa T, Konnai S, Nishimori A et al. Bovine immunoinhibitory receptors contribute to the suppression of Mycobacterium avium subsp.paratuberculosis-specific T-cell responses Infect Immun. 2015-10-19 [PMID: 26483406] (FLOW)

Freeman CM, Crudgington S et al. Design of a multi-center immunophenotyping analysis of peripheral blood, sputum and bronchoalveolar lavage fluid in the Subpopulations and Intermediate Outcome Measures in COPD Study (SPIROMICS). J Transl Med 2015-01-01 [PMID: 25622723] (FLOW)

Robinson AP, Rodgers JM, Goings GE, Miller SD. Characterization of Oligodendroglial Populations in Mouse Demyelinating Disease Using Flow Cytometry: Clues for MS Pathogenesis. PLoS One 2014-01-01 [PMID: 25247590] (FLOW)

Taichman RS, Patel LR, Bedenis R et al. GAS6 Receptor Status Is Associated with Dormancy and Bone Metastatic Tumor Formation. PLoS One 2013-01-01 [PMID: 23637920] (FLOW)

Jones CL, Ferreira S, McKenzie RC et al. Regulation of T-Plastin Expression by Promoter Hypomethylation in Primary Cutaneous T-Cell Lymphoma. J Invest Dermatol 2013-01-01 [PMID: 22495182] (FLOW)

Kim N, Kukkonen S, Gupta S, Aldovini A. Association of Tat with Promoters of PTEN and PP2A Subunits Is Key to Transcriptional Activation of Apoptotic Pathways in HIV-Infected CD4+ T Cells. PLoS Pathog. 2010-01-01 [PMID: 20862322] (FLOW)



### **Novus Biologicals USA**

10730 E. Briarwood Avenue  
Centennial, CO 80112

USA

Phone: 303.730.1950

Toll Free: 1.888.506.6887

Fax: 303.730.1966

[nb-customerservice@bio-techne.com](mailto:nb-customerservice@bio-techne.com)

### **Bio-Techne Canada**

21 Canmotor Ave

Toronto, ON M8Z 4E6

Canada

Phone: 905.827.6400

Toll Free: 855.668.8722

Fax: 905.827.6402

[canada.inquires@bio-techne.com](mailto:canada.inquires@bio-techne.com)

### **Bio-Techne Ltd**

19 Barton Lane

Abingdon Science Park

Abingdon, OX14 3NB, United Kingdom

Phone: (44) (0) 1235 529449

Free Phone: 0800 37 34 15

Fax: (44) (0) 1235 533420

[info.EMEA@bio-techne.com](mailto:info.EMEA@bio-techne.com)

### **General Contact Information**

[www.novusbio.com](http://www.novusbio.com)

Technical Support: [nb-technical@bio-techne.com](mailto:nb-technical@bio-techne.com)

Orders: [nb-customerservice@bio-techne.com](mailto:nb-customerservice@bio-techne.com)

General: [novus@novusbio.com](mailto:novus@novusbio.com)

### **Limitations**

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Kits are guaranteed for 6 months from date of receipt.

For more information on our 100% guarantee, please visit [www.novusbio.com/guarantee](http://www.novusbio.com/guarantee)

Earn gift cards/discounts by submitting a review: [www.novusbio.com/reviews/submit/765-0015](http://www.novusbio.com/reviews/submit/765-0015)

Earn gift cards/discounts by submitting a publication using this product:  
[www.novusbio.com/publications](http://www.novusbio.com/publications)

