

Product Datasheet

Lightning-Link (R) PE-Atto 594 Antibody Labeling Kit 768-0030

Unit Size: 3 x 10ug Reaction

Store at -20C.

www.novusbio.com



technical@novusbio.com

Publications: 1

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/768-0030

Updated 8/19/2024 v.20.1

Earn rewards for product
reviews and publications.

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/768-0030



768-0030**Lightning-Link (R) PE-Atto 594 Antibody Labeling Kit**

Product Information	
Unit Size	3 x 10ug Reaction
Concentration	Concentration is not relevant for this product. Please see the protocols for proper use of this product.
Storage	Store at -20C.
Conjugate	PE/Atto594
Product Description	
Description	<p>Lightning-Link antibody labeling kits enable the direct labeling of antibodies, proteins, peptides or other biomolecules for use in R&D applications, drug discovery and the development of diagnostic kits (See protocol for further information).</p> <p>Our PE/Atto594 antibody labeling kit enables the direct conjugation of the PE/Atto594 tandem label to any biomolecule with an available amine group. The researcher simply pipettes the antibody or other biomolecule into the vial of Lightning-Link label and incubates for 3 hours.</p> <p>FeaturesQuick and easy to use BenefitsSave time, no special knowledge required No separation steps 100% recovery - no antibody/protein loss Can be used in a wide range of applications Flexible Freeze dried Ships at ambient temperature, long shelf-life Fully scalable (10 ug to 1 g or more) Easy transfer from R&D to manufacturing Stringently QC tested Consistent high quality, excellent batch-to-batch reproducibility Large number of labels available Experimental flexibility Reliable: nearly 300 references Successfully used in many fields of research</p> <p>The PE has three maximal absorbance values of 498, 544 and 566nm with the optimal depending on the application. The PE functions as an energy donor for the Atto594 fluorescent dye and transfers its energy via energy resonance transfer. The Atto594 emits the energy received from the PE in the form of long wavelength light at 632nm</p> <p>Learn more about Lightning-Link™ Conjugation Kits by reading FAQs</p> <p>For more information please check out these useful links! Antibody Labeling Guide Antibody Conjugation Illustrated Assay</p>
Kit Components	1 or 3 glass vial(s) of Lightning-Link mix, 1 vial of LL-Modifier reagent, 1 vial of LL-Quencher reagent
Notes	<p>This product is manufactured by Abcam and distributed by Novus Biologicals.</p> <p>This product is for research use only and is not approved for use in humans or in clinical diagnosis. This product is guaranteed for 1 year from date of receipt and this statement overrides any mentioned guarantee period on the limitations section of this products datasheet. Please contact technical@novusbio.com with questions.</p>
Product Application Details	

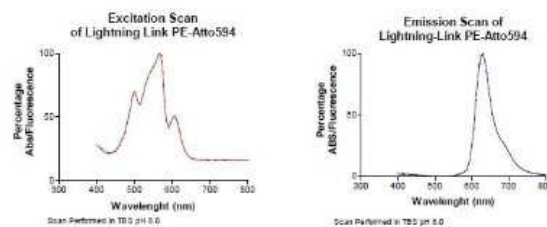


Application Notes

By circumventing the desalting or dialysis steps that commonly interrupt traditional antibody conjugation procedures, LightningLink technology can be used to label both small (e.g. 10 ug) and large quantities of primary antibodies with ease. Batch-to-batch variation upon scale up is minimal as the process is so simple, and recoveries are always 100%.

Images

Lightning-Link PE-Atto 594 Antibody Labeling Kit [768-0030]

**Publications**

Luk HY, McFarlin BK, Vingren JL. Using image-based flow cytometry to monitor satellite cells proliferation and differentiation in vitro. Methods. 2017-01-01 [PMID: 27524556]





Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112
USA
Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave
Toronto, ON M8Z 4E6
Canada
Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB, United Kingdom
Phone: (44) (0) 1235 529449
Free Phone: 0800 37 34 15
Fax: (44) (0) 1235 533420
info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com
Technical Support: nb-technical@bio-techne.com
Orders: nb-customerservice@bio-techne.com
General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Kits are guaranteed for 6 months from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/768-0030

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

