## **Product Datasheet**

# Aurora B Antibody - BSA Free NB100-294

Unit Size: 0.1 ml

Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.

www.novusbio.com



technical@novusbio.com

Reviews: 5 Publications: 12

Protocols, Publications, Related Products, Reviews, Research Tools and Images at: www.novusbio.com/NB100-294

Updated 10/23/2024 v.20.1

Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications Submit a review at www.novusbio.com/reviews/destination/NB100-294



## NB100-294

Aurora B Antibody - BSA Free


| Aurora B Antibody - BSA Free |   |
|------------------------------|---|
| Product Information          |   |
| Unit Size                    | 0.1 ml  |
| Concentration                | 12.5 mg/ml  |
| Storage                      | Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.  |
| Clonality                    | Polyclonal  |
| Preservative                 | 0.09% Sodium Azide  |
| Isotype                      | IgG   |
| Purity                       | IgG purified  |
| Buffer                       | 10mM PBS (pH 7.4)   |
| Target Molecular Weight      | 41 kDa  |
| Product Description          |   |
| Host                         | Rabbit  |
| Gene ID                      | 9212  |
| Gene Symbol                  | AURKB   |
| Species                      | Human, Mouse, Rat   |
| Marker                       | Mitosis Marker  |
| Specificity/Sensitivity      | Aurora B  |
| Immunogen                    | Synthetic peptide corresponding to amino acids 1-19 of human Aurora B with C-terminal added cysteine conjugated to KLH. Uniprot accession no. Q96GD4  |
| Product Application Details  |   |
| Applications                 | Western Blot, Immunocytochemistry/ Immunofluorescence, Immunohistochemistry, Immunohistochemistry-Paraffin, Immunoprecipitation, Microarray   |
| Recommended Dilutions        | Western Blot 1:1000-1:2000, Immunohistochemistry 1:10-1:500, Immunocytochemistry/ Immunofluorescence 1:50, Immunoprecipitation 50-60 ug, Immunohistochemistry-Paraffin 1:10-1:500, Microarray |
| Application Notes            | By Western blot, it detects a band of approximately 41 kDa. Detection of the Aurora B band by immunoblotting is specifically inhibited with the immunizing                                    |




peptide.

## **Images**

Western Blot: Aurora B Antibody - BSA Free [NB100-294] - Whole extract of PC-12 cells were separated on SDS-PAGE and probed with Anti-Aurora B antibody produced in Rabbit. The antibody was developed using 1:10,000 Anti-Rabbit IgG (whole molecule)-Peroxidase antibody produced in Goat. Lanes: 1. 1:1000 antibody 2. 1:2000 antibody 3. Negative Control

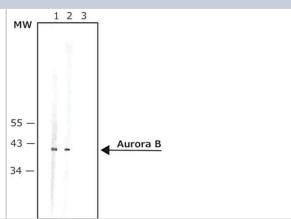

Immunocytochemistry/Immunofluorescence: Aurora B Antibody - BSA Free [NB100-294] - HeLa human cervix epitheloid carcinoma cells were fixed and permeabilized with 4% paraformaldehyde followed by cold methanol. Fixed cells were stained with 1:50 Anti-Aurora B antibody produced in rabbit. The antibody was developed using Anti-Rabbit IgG (whole molecule)-FITC antibody produced in goat.

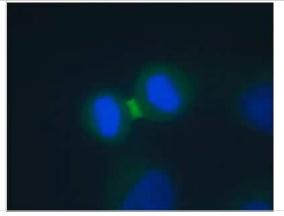

### **Publications**

Carlotta Penzo, Laura Arnoldo, Silvia Pegoraro, Sara Petrosino, Gloria Ros, Rossella Zanin, Jacek R. Wiśniewski, Guidalberto Manfioletti, Riccardo Sgarra HMGA1 Modulates Gene Transcription Sustaining a Tumor Signalling Pathway Acting on the Epigenetic Status of Triple-Negative Breast Cancer Cells Cancers 2019-08-02 [PMID: 31382504]

Chang X Hi-511 Inhibits Both BRAF and AURKB to Overcome Malignant Melanoma Drug Resistance and Metastasis Thesis 2020-01-01

Chang X, Zhang T, Wang Q et al. HI-511 overcomes melanoma drug resistance via targeting AURKB and BRAF V600E Theranostics 2020-08-01 [PMID: 32863956] (Human)

Fang H, Niu K, Mo D et al. RecQL4-Aurora B kinase axis is essential for cellular proliferation, cell cycle progression, and mitotic integrity. Oncogenesis 2018-09-12 [PMID: 30206236] (ICC/IF, Human)

Dumitru AMG, Rusin SF, Clark AEM et al. Cyclin A/Cdk1 modulates Plk1 activity in prometaphase to regulate kinetochore-microtubule attachment stability Elife. 2017-11-20 [PMID: 29154753] (ICC/IF, WB, Human)

#### Details:

HEC1 IF antibody used at a dilution of 1:500

Sijare F, Geibler AL, Fichter CD et al. Aurora B expression and histone variant H1.4S27 phosphorylation are no longer coordinated during metaphase in aneuploid colorectal carcinomas Virchows Arch. 2015-02-14 [PMID: 25680570] (ICC/IF)

Kleyman M, Kabeche L, Compton DA. STAG2 promotes error correction in mitosis by regulating kinetochore-microtubule attachments. J. Cell. Sci. 2014-07-29 [PMID: 25074805] (ICC/IF, Human)

### Details:

Outer kinetochore component Hec1/NDC80 antibody used in ICC-IF on HCT116 and RPE1 cells for measuring centromere stretch (Fig. 1D, supplementary material Fig. S2A)

Kabeche L, Compton DA. Cyclin A regulates kinetochore microtubules to promote faithful chromosome segregation. Nature. 2013-10-03 [PMID: 24013174] (ICC/IF, Human)

Kabeche L, Compton DA. Checkpoint-independent stabilization of kinetochore-microtubule attachments by Mad2 in human cells. Curr Biol 2012-04-01 [PMID: 22405866] (ICC/IF, Human)

Lee EC, Frolov A, Li R et al. Targeting Aurora kinases for the treatment of prostate cancer. Cancer Res 2006-05-01 [PMID: 16707419]

Beussel S, Hasenburg A, Bogatyreva L, Hauschke D, Werner M, Lassmann S. Aurora-B protein expression is linked to initial response to taxane-based first-line chemotherapy in stage III ovarian carcinoma. J Clin Pathol;65(1):29-35. 2012-01-01 [PMID: 22011448] (IF/IHC, Human)

Bufo P, Sanguedolce F, Tortorella S et al. Expression of mitotic kinases phospho-aurora A and aurora B correlates with clinical and pathological parameters in bladder neoplasms. Histol Histopathol;25(11):1371-7. 2010-11-01 [PMID: 20865660]





## Novus Biologicals USA

10730 E. Briarwood Avenue Centennial, CO 80112

USA

Phone: 303.730.1950 Toll Free: 1.888.506.6887

Fax: 303.730.1966

nb-customerservice@bio-techne.com

## **Bio-Techne Canada**

21 Canmotor Ave Toronto, ON M8Z 4E6

Canada

Phone: 905.827.6400 Toll Free: 855.668.8722 Fax: 905.827.6402

canada.inquires@bio-techne.com

## **Bio-Techne Ltd**

19 Barton Lane Abingdon Science Park Abingdon, OX14 3NB, United Kingdom Phone: (44) (0) 1235 529449

Free Phone: 0800 37 34 15 Fax: (44) (0) 1235 533420 info.EMEA@bio-techne.com

## **General Contact Information**

www.novusbio.com

Technical Support: nb-technical@bio-

techne.com

Orders: nb-customerservice@bio-techne.com

General: novus@novusbio.com

## **Products Related to NB100-294**

NBL1-07860 Aurora B Overexpression Lysate

HAF008 Goat anti-Rabbit IgG Secondary Antibody [HRP]

NB7160 Goat anti-Rabbit IgG (H+L) Secondary Antibody [HRP]

NBP2-24891 Rabbit IgG Isotype Control

#### Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NB100-294

Earn gift cards/discounts by submitting a publication using this product: www.novusbio.com/publications

