

Product Datasheet

Claudin-5 Antibody NB100-91713

Unit Size: 0.1 ml

Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.

www.novusbio.com



technical@novusbio.com

Publications: 2

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/NB100-91713

Updated 1/25/2018 v.20.1

Earn rewards for product
reviews and publications.

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/NB100-91713



NB100-91713

Claudin-5 Antibody


Product Information	
Unit Size	0.1 ml
Concentration	1.0 mg/ml
Storage	Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.
Clonality	Polyclonal
Preservative	0.05% Sodium Azide
Purity	Immunogen affinity purified
Buffer	PBS (pH 7.2)
Target Molecular Weight	23 kDa
Product Description	
Host	Rabbit
Gene ID	7122
Gene Symbol	CLDN5
Species	Human, Mouse, Rat
Reactivity Notes	Human, Mouse and Rat.
Immunogen	Synthetic peptide, corresponding to amino acids 171-220 of Human Claudin-5.
Product Application Details	
Applications	Western Blot, Immunocytochemistry/Immunofluorescence, Immunohistochemistry, Immunohistochemistry-Paraffin
Recommended Dilutions	Western Blot 1:500-1:1000, Immunohistochemistry 1:10-1:500, Immunocytochemistry/Immunofluorescence 1:50-1:200, Immunohistochemistry-Paraffin 1:50-1:200
Application Notes	Western blot (WB) analysis of Claudin 5 (S201) pAb in extracts from A549 cells. Immunohistochemistry (IHC) analyzes of Claudin 5 (S201) pAb in paraffin-embedded human brain tissue.

Images

Western Blot: Claudin 5 Antibody [NB100-91713] - Extracts from A549 cells



Immunohistochemistry-Paraffin: Claudin 5 Antibody [NB100-91713] - Paraffin-embedded human brain tissue.

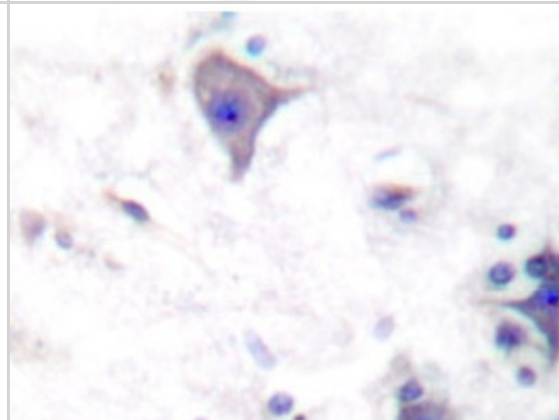

Publications

Ding R, Chen Y, Yang S et al. Blood-brain barrier disruption induced by hemoglobin in vivo: Involvement of up-regulation of nitric oxide synthase and peroxynitrite formation. *Brain Res.* 2014 May 09 [PMID: 24814387] (IHC-P, WB, Rat)

Madlambayan GJ, Butler JM, Hosaka K et al. Bone marrow stem progenitor cell contribution to neovasculogenesis is dependent on model system with SDF-1 as a permissive trigger. *Blood* :blood-2009-2003-211342. 2009 [PMID: 19717647]

Procedures

Western Blot Protocol for Claudin 5 Antibody (NB100-91713)

Western Blot Protocol:

1. Prepare cell / tissue lysate(s) (see cell / tissue lysis protocol)
2. Make up lysates to a total volume of 10-20 ul with 4x SDS sample buffer (see recipes) containing a reducing agent (either 10% B-ME or 0.3M DTT in the 4x SB).
3. Prepare SDS-PAGE gel (see SDS-PAGE gel pouring).
4. Place gel in running chamber and fill with 1 x SDS Running Buffer. Use 10x stock (see recipes).
5. Heat samples at 95-100C for 5 minutes.
6. Load 5ul of a control pre-stained protein ladder to track sample migration in the first lane of the gel. This is also useful for the transfer step to verify the procedure is successful. These are widely available for most suppliers (Invitrogen, Biorad, Fermentas, etc).
7. Load boiled samples. Electrophorese at constant 100V until the blue bromophenol dye reach the bottom of the gel.

Protein Transfer:

1. Place gel and blotting pads in transfer buffer solution.
2. Transfer gel to filter paper wetted in transfer buffer solution.
3. Assemble transfer sandwich by orientating cathode, filter paper, gel, membrane (nitrocellulose / PVDF), filter paper, anode so protein transfer goes in the direction of cathode to anode. These are typically coloured black (cathode) and red (anode). The choice of membrane (nitrocellulose / PVDF) depends on the protein studied and the antibodies used. Generally you will get the use instructions you need through antibody guidelines. PVDF is particularly useful if you are working with small or highly charged proteins. Remember to make PVDF soaked in methanol for 5 minutes prior to use.
4. Electrophorese for 90 minutes at 100V or overnight at 4 degrees C at 20V. Proceed electrophoresis on ice or in a cold room as the transfer generates a lot of heat!
5. Check membrane for transfer of control protein ladder. If using nitrocellulose, sample transfer and relative concentration can be assessed by staining the membrane with Ponceau S solution (0.1% ponceau, 5% acetic acid). This is a reversible stain which is removed by a simple wash in TBST.

Protein Labeling:

1. Incubate membrane in blocking solution for 1 hour at room temperature or at 4 degrees C overnight with shaking. The blocking solution is normally composed of 5% non-fat milk in TBS-T, and some antibodies require BSA in place of milk. This is normally clear in the manufacturers instructions for the antibody for testing. To make TBS-T, using 10 x TBS, add Tween-20 to a final concentration of 0.05%.
2. Incubate primary antibody overnight or at room temperature for 2 hours prepared in 5% milk/TBST.
3. Remove antibody solution and wash membrane three times for 5-15 minutes in TBST.
4. Incubate membrane with an appropriate secondary antibody (peroxidase conjugated) for 1 hour at room temperature in 5% milk/TBST. For example if your primary antibody was raised in rabbit, a goat anti-rabbit HRP should be used as a secondary antibody.
5. Remove antibody solution and wash membrane three times for 5-15 minutes in TBST.

Protein Detection:

1. There are numerous commercial available chemiluminescence reagents (Amersham, Pierce, Invitrogen) and a range of sensitivities of detection levels is included. These typically take the form of two solutions which are combined and then incubated immediately on the membrane for 1 - 5 minutes.
2. Make membrane exposed to X-ray film for 1 minute to 1 hour, depending on protein signal and chemiluminescence method.



Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112
USA

Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
novus@novusbio.com

Novus Biologicals Canada

461 North Service Road West, Unit B37
Oakville, ON L6M 2V5
Canada

Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada@novusbio.com

Novus Biologicals Europe

19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB, United Kingdom

Phone: (44) (0) 1235 529449
Free Phone: 0800 37 34 15
Fax: (44) (0) 1235 533420
info@bio-techne.com

General Contact Information

www.novusbio.com
Technical Support: technical@novusbio.com
Orders: orders@novusbio.com
General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NB100-91713

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

