

Product Datasheet

beta-Arrestin 1 Antibody (E274)

NB110-55485

Unit Size: 0.1 ml

Store at -20C. Avoid freeze-thaw cycles.

www.novusbio.com



technical@novusbio.com

Publications: 3

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/NB110-55485

Updated 10/3/2017 v.20.1

Earn rewards for product
reviews and publications.

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/NB110-55485



NB110-55485

beta-Arrestin 1 Antibody (E274)

Product Information

| | |
|--------------------------------|---|
| Unit Size | 0.1 ml |
| Concentration | Please see the vial label for concentration. If unlisted please contact technical services. |
| Storage | Store at -20C. Avoid freeze-thaw cycles. |
| Clonality | Monoclonal |
| Clone | E274 |
| Preservative | 0.01% Sodium Azide |
| Isotype | IgG |
| Purity | Protein A or G purified |
| Buffer | 49% PBS, 0.05% BSA and 50% Glycerol |
| Target Molecular Weight | 47 kDa |

Product Description


| | |
|-------------------------|---|
| Host | Rabbit |
| Gene ID | 408 |
| Gene Symbol | ARRB1 |
| Species | Human, Mouse, Rat |
| Reactivity Notes | Predicted to work with: Bovine |
| Immunogen | A synthetic peptide corresponding to residues near N-terminus of human - arrestin-1 was used as immunogen. |
| Notes | Produced using Abcam's RabMab® technology. RabMab® technology is covered by the following U.S. Patents, No. 5,675,063 and/or 7,429,487. |

Product Application Details


| | |
|------------------------------|--|
| Applications | Western Blot, Flow Cytometry, Immunocytochemistry/Immunofluorescence, Immunohistochemistry, Immunohistochemistry-Paraffin |
| Recommended Dilutions | Western Blot 1:1000-10000, Flow Cytometry 1:40, Immunohistochemistry 1:10-1:500, Immunocytochemistry/Immunofluorescence 1:100, Immunohistochemistry-Paraffin 1:100 |

Images


Western Blot: beta-Arrestin 1 Antibody (E274) [NB110-55485] - 293 cell lysate.

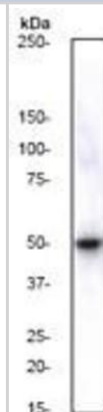


Immunocytochemistry/Immunofluorescence: beta-Arrestin 1 Antibody (E274) [NB110-55485] - Human prostate cancer cell line.

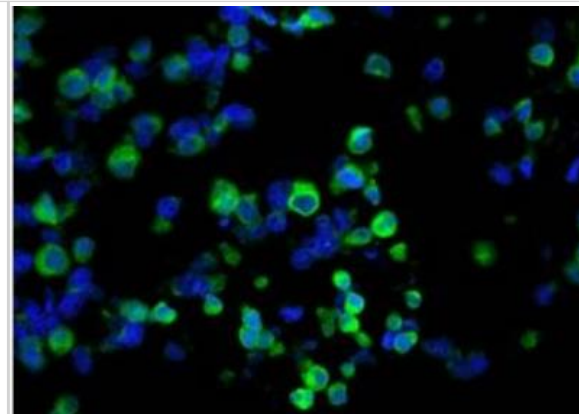


Immunohistochemistry-Paraffin: beta-Arrestin 1 Antibody (E274) [NB110-55485] - Human prostate carcinoma tissue section.

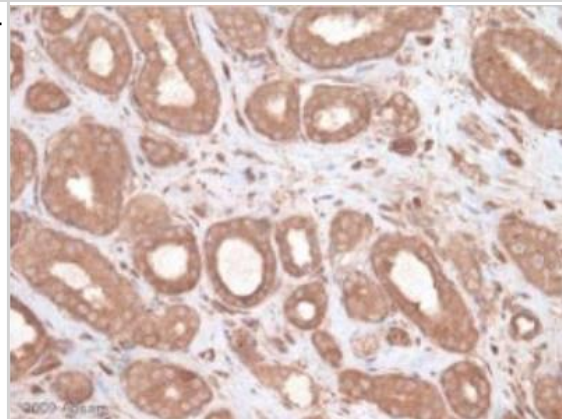

Flow Cytometry: beta-Arrestin 1 Antibody (E274) [NB110-55485] - Overlay histogram showing PC3 cells stained with beta Arrestin 1 antibody [E274] (red line).

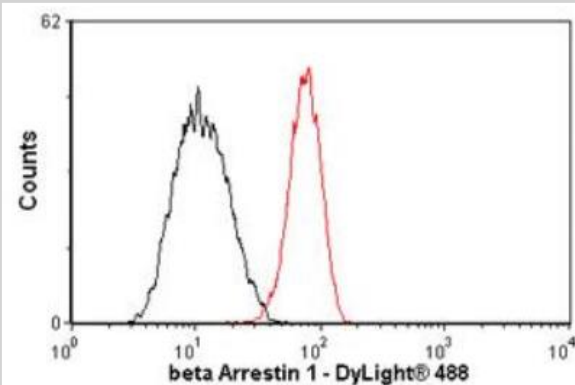

Immunohistochemistry: beta-Arrestin 1 Antibody (E274) [NB110-55485] - Immunohistochemical analysis of paraffin-embedded human lymph node using anti-β-arrestin

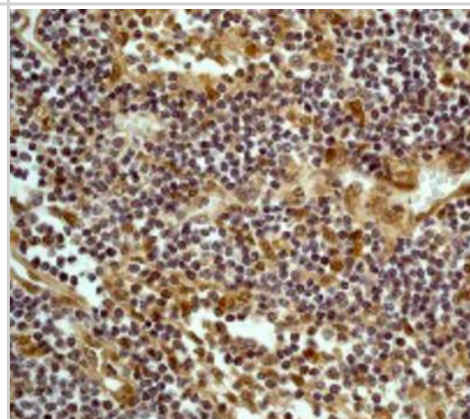

Publications

Kumar R, Samal SK, Routray S et al. Identification of oral cancer related candidate genes by integrating protein-protein interactions, gene ontology, pathway analysis and immunohistochemistry. Sci Rep May 30 2017 12:00AM [PMID: 28559546] (IHC, Human)

Wang WC, Schillinger RM, Malone MM, Liggett S. Paradoxical attenuation of B2-AR function in airway smooth muscle by Gi-mediated counterregulation in transgenic mice overexpressing type 5 adenylyl cyclase Am J Physiol Lung Cell Mol Physiol 2011 Mar [PMID: 21131397] (WB, Mouse)

Oakley RH, Revollo J, Cidlowski JA. Glucocorticoids regulate arrestin gene expression and redirect the signaling profile of G protein-coupled receptors Proc Natl Acad Sci U S A 2012 Oct 8 [PMID: 23045642] (WB, ICC/IF, Human)



Procedures



Immunohistochemistry Protocol for Arrestin beta 1 Antibody (NB110-55485)**Immunohistochemistry Protocol for Paraffin-embedded Tissues****1. Solutions and reagents****1.1. Xylene****1.2. Ethanol, anhydrous denatured, histological grade (100%, 95%, 70%)****1.3. Washing buffer:**

TBST washing buffer: 1XTBS/0.1% Tween-20

To prepare stock solution of 10X TBS: add 24.2 g Trizma base and 80 g sodium chloride to 1L of dH₂O. Adjust pH to 7.6.

Working solution. 1XTBST/0.1% Tween-20: add 100ml 10XTBS to 900 ml dH₂O. Add 1 ml Tween-20 and mix well.

1.4. Distilled water (dH₂O)**1.5. Antigen Retrieval Solution:**

0.01M Sodium Citrate Buffer, pH 6.0

To prepare stock solutions:

Solution A. 0.1 M citric acid solution: dissolve 21.0 g of citric acid, monohydrate (C₆H₈O₇·H₂O) in 1 liter of dH₂O.

Solution B. 0.1M sodium citrate solution: dissolve 29.4 g trisodium citrate dihydrate (C₆H₅Na₃O₇·2H₂O) in 1 liter of dH₂O.

Working solution: Add 9 ml of Stock solution A and 41 ml of stock solution B to 450 ml of dH₂O. Adjust pH to 6.0.

1.6. 3% Hydrogene Peroxide**1.7. Blocking buffer:**

PBS (Dulbeccos Phosphate Buffered Salts, 1X, catalog #21-031-CV from Mediatech, Inc.) + 10% serum (serum origin depends on the host of the secondary antibody)

1.8. Hematoxylin QS (catalog #H-3404 from Vector Laboratories, Inc.)**1.9. Permanent Mounting medium (VectaMount, catalog# H-5000 Vector Laboratories, Inc.)****2. Protocol****2.1. Deparaffinization/Rehydration****2.1.1. Heat slides in an oven at 65C for 1 hour.**

2.1.2. De-paraffinize/hydrate using the following series of washes: two Xylene washes (5 min each), followed by two 100% ethanol rinses (5 min each), followed by 95% ethanol, 70% ethanol, 50% ethanol, 30% ethanol, followed by H₂O and a TBST wash for 5 min on a shaker.

2.2. Antigen Retrieval

2.2.1. Immerse slides into staining dish containing Antigen Retrieval Solution.

2.2.2. Place covered staining dish into the rice cooker. Add 120 ml of dH₂O.

2.2.3. When cook is turned to warm (about 20 to 30 min), unplug the cooker and remove the staining dish to the bench top.

2.2.4. Allow to cool down, without cover, for 20 min.

2.3. Staining

2.3.1. Wash slides with TBST for 5 min on a shaker.

2.3.2. Inactivate endogenous peroxidase by covering tissue with 3% hydrogen peroxide for 10 min.

2.3.3. Wash slides three times with TBST (3 min each on a shaker).

2.3.4. Block slides with the blocking solution for 1 hour.

2.3.5. Dilute primary antibody in the blocking buffer per recommendation on the data sheet.

2.3.6. Apply primary antibody to each section and incubate overnight in the humidified chamber (4C).

2.3.7. Wash slides three times with TBST (3 min each on a shaker).

2.3.8. Apply to each section secondary HRP-conjugated anti-rabbit antibody diluted in the blocking solution per manufacturers recommendation; incubate for 1 hour at room temperature.

2.3.9. Wash slides three times with TBST (3 min each on a shaker).

2.3.10. Add freshly prepared DAB substrate to the sections.

2.3.11. Incubate tissue sections with the substrate at room temperature until suitable staining develops (generally 2 to 5 min).

2.3.12. Rinse sections with water.

2.3.13. Counterstain with Hematoxylin.

2.3.14. Rinse sections with water.

2.3.15. Dehydrate samples using two rinses with 100% Ethanol (20 dips per rinse) followed by two rinses with Xylene (30 dips per rinse).

2.3.16. Mount coverslips on slides using Permount medium.



Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112
USA
Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
novus@novusbio.com

Novus Biologicals Canada

461 North Service Road West, Unit B37
Oakville, ON L6M 2V5
Canada
Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada@novusbio.com

Novus Biologicals Europe

19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB, United Kingdom
Phone: (44) (0) 1235 529449
Free Phone: 0800 37 34 15
Fax: (44) (0) 1235 533420
info@bio-techne.com

General Contact Information

www.novusbio.com
Technical Support: technical@novusbio.com
Orders: orders@novusbio.com
General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NB110-55485

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

