Product Datasheet

Coronin-1a Antibody - BSA Free NB110-58367

Unit Size: 0.5 ml

Store at 4C in the dark.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at: www.novusbio.com/NB110-58367

Updated 2/21/2025 v.20.1

Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications
Submit a review at www.novusbio.com/reviews/destination/NB110-58367



NB110-58367


Coronin-1a Antibody - BSA Free

Coronin-1a Antibody - BSA Free	
Product Information	
Unit Size	0.5 ml
Concentration	0.1 mg/ml
Storage	Store at 4C in the dark.
Clonality	Polyclonal
Preservative	0.02% Sodium Azide
Isotype	IgY
Purity	Immunogen affinity purified
Buffer	10mM PBS (0.9% isotonic, w/v, pH 7.2)
Target Molecular Weight	51 kDa
Product Description	
Host	Chicken
Gene ID	11151
Gene Symbol	CORO1A
Species	Human, Mouse, Rat
Marker	Microglial Cell Marker
Immunogen	Chickens were immunized with a synthetic peptide/keyhole limpet hemocyanin (KLH) conjugate. This synthetic peptide corresponded to a sequence shared between the human (NM_007074) and mouse (NM_009898) gene products.
Notes	Chicken products cannot be exported to Canada. Purification Notes After repeated injections, immune eggs were collected, the IgY fractions were purified from the yolks, and then affinity purified using a peptide column. The concentration of the eluate was then adjusted to 100 ug/ml, and the preparation was filter-sterilized. Storage Notes Store at 4C in the dark. Under these conditions, the antibodies should have a shelf life of at least 12 months (provided they remain sterile). Do not freeze these antibodies unless you want to store them for longer periods of time. Note, however, that each time an antibody preparation is frozen, about half of its binding activity is lost.
Product Application Details	
Applications	Western Blot, Immunocytochemistry/ Immunofluorescence, Immunohistochemistry
Recommended Dilutions	Western Blot 1:5000-1:10000, Immunohistochemistry 1:1000-1:2000, Immunocytochemistry/ Immunofluorescence 1:1000-1:2000



Images

Immunocytochemistry/Immunofluorescence: Coronin-1a Antibody [NB110-58367] - Anti-Coronin 1A

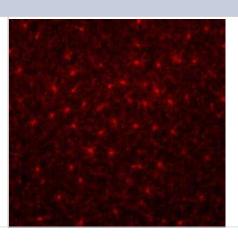

Procedures

Immunohistochemistry Chicken IgY Protocol (NB110-58367)

Immunohistochemistry Chicken IgY Protocol (NB110-58367): Citrate Buffer Antigen Retrieval Protocol

necessary to compensate for their larger size.

Background: Formaldehyde fixation (2% or 4%, or as a component of 10% formalin) produces protein cross-links in tissues that tends to interfere with antibody penetration. This seems to be particularly true of paraffin- embedded formaldehyde-fixed tissue. Since chicken IgY antibodies are larger than rabbit or mouse IgG's, "extra steps" may be

The citrate-based "antigen retrieval" protocol outlined below has been shown to improve chicken IgY antibody penetration into 4% formalde- hyde-fixed paraffin-embedded sections, and can increase the degree and intensity of immunoreactivity and immunostaining.

Reagents (NOTE: You can use either the Sodium Citrate or Citric Acid Buffers in step #3, below)

"Sodium Citrate Buffer" (10mM Sodium Citrate, 0.05% Tween 20, pH 6.0)

Weigh out 2.94 grams of trisodium citrate (dihydrate). Dissolve in approximately 900 mls of deionized, distilled water. Adjust the pH to 6.00 with 1.0 N HCl. Add 0.5 ml of Tween-20. Mix. Bring up the volume to 1.0 litres with water. Store this solution at room temperature for 3 months or at 4C for longer periods.

"Citric Acid Buffer" (10mM Citric Acid, 0.05% Tween 20, pH 6.0)

Weigh out 1.92 grams of citric acid (anhydrous). Dissolve in approximately 900 mls of deionized, distilled water. Adjust the pH to 6.0 with 1.0 N NaOH. Add 0.5 ml of Tween-20. Mix. Bring up the volume to 1.0 litres with water. Store this solution at room temperature for 3 months or at 4C for longer periods.

"Phosphate-Buffered Saline" [PBS, 10 mM Sodium phosphate-buffered (pH 7.2) isotonic (0.9%, w/v) saline solution] PBS Tween (0.05% Tween 20 in PBS) Ethanol (80%, 90%, 95%, 100%) diluted with water

Xylene

Procedure (for use with paraffin-embedded sections):

- 1 Deparaffinize tissue sections in 2 changes of xylene (5 minutes each).
- 2. Hydrate in 2 changes of 100% ethanol (3 minutes each), 95% ethanol (1 minute), 90% ethanol (1 minute), 80% ethanol (1 minute). Rinse in distilled water.
- 3. Pre-heat steamer or water bath with staining dish containing either Sodium Citrate Buffer or Citrate Buffer. Wait until temperature reaches 95-100 degrees C.

NOTE: Microwave or pressure cooker can be used as an alternative as a heating source.

- 4. Immerse slides in the staining dish. Place the lid loosely on the staining dish and incubate for 20-40 minutes (optimal incubation times will vary).
- 5. Remove the staining dish, and allow it to cool to room temperature (for 20 minutes or so).



Rinse sections in PBS Tween twice for 2 minutes each time.

NOTE: The remainder of this protocol is meant to be a suggestion, and can be substituted with your regular immunostaining protocol.

- 7. Block sections for 30 minutes with Blocking buffer diluted 1:10 with water.
- 8. Incubate sections with primary antibody at appropriate dilution in antibody dilution buffer overnight at 4 degrees C. Since chicken IgY antibodies are larger than mammalian IgG's, this overnight incubation allows more time for antibody penetration into tissue sections.
- 9. Rinse sections with PBS Tween 20 twice for 5 minutes each time.
- 10. Incubate sections with labeled secondary antibody (see NOTE, below) at appropriate dilution (for one hour at room temperature) in a 1:100 dilution of blocking buffer (diluted in PBS).
- Rinse with PBS Tween 20 for three times for 5 minutes each time.

NOTE: This protocol may use HRP- or fluorescently-labeled secondary antibodies produced in goats or rabbits.

References:

- 1. Shi SR, Chaiwun B, Young L, Cote RJ, Taylor CR. (1993). Antigen retrieval technique utilizing citrate buffer or urea solution for immunohistochemical demonstration of androgen receptor in formalin-fixed paraffin sections. J Histochem Cytochem 41 (11): 1599-1604.
- 2. Kanai K, Nunoya T, Shibuya K, Nakamura T, Tajima M (1998). Variations in effectiveness of antigen retrieval pretreatments for diagnostic immunohistochemistry. Res Vet Sci 64 (1): 57-61.
- 3. Brown RW, Chirala R. (1995). Utility of microwave-citrate antigen retrieval in diagnostic immunohistochemistry. Mod Pathol 8 (5): 515-20.
- 4. Morgan JM, Navabi H, Schmid KW, Jasani B (1994). Possible role of tissue-bound calcium ions in citrate-mediated high-temperature antigen retrieval. J Pathol 174 (4): 301-7.
- 5. Pellicer EM, Sundblad A (1994). Antigen retrieval by microwave oven with buffer of citric acid. Medicina (B Aires). 54 (2): 129-32.
- 6. Shi SR, Chaiwun B, Young L, Cote RJ, Taylor CR (1993). Antigen retrieval technique utilizing citrate buffer or urea solution for immunohistochemical demonstration of androgen receptor in formalin-fixed paraffin sections. J Histochem Cytochem 41 (11): 1599-604.





Novus Biologicals USA

10730 E. Briarwood Avenue Centennial, CO 80112

USA

Phone: 303.730.1950 Toll Free: 1.888.506.6887

Fax: 303.730.1966

nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave Toronto, ON M8Z 4E6

Canada

Phone: 905.827.6400 Toll Free: 855.668.8722 Fax: 905.827.6402

canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane Abingdon Science Park Abingdon, OX14 3NB, United Kingdom Phone: (44) (0) 1235 529449

Free Phone: 0800 37 34 15 Fax: (44) (0) 1235 533420 info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com

Technical Support: nb-technical@bio-

techne.com

Orders: nb-customerservice@bio-techne.com

General: novus@novusbio.com

Products Related to NB110-58367

BAF010 Goat anti-Chicken IgY Secondary Antibody [Biotin]

NB7276 Goat anti-Chicken IgM Heavy Chain Secondary Antibody
H00011151-Q01-10ug Recombinant Human Coronin-1a GST (N-Term) Protein

NBL1-09395 Coronin-1a Overexpression Lysate

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NB110-58367

Earn gift cards/discounts by submitting a publication using this product: www.novusbio.com/publications

