

Product Datasheet

TRAP1 Antibody (TRAP1-6)

NB300-555

Unit Size: 100ug

Store at -20C. Avoid freeze-thaw cycles.

www.novusbio.com



technical@novusbio.com

Publications: 4

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/NB300-555

Updated 10/23/2024 v.20.1

Earn rewards for product
reviews and publications.

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/NB300-555



NB300-555**TRAP1 Antibody (TRAP1-6)**


Product Information	
Unit Size	100ug
Concentration	1 mg/ml
Storage	Store at -20C. Avoid freeze-thaw cycles.
Clonality	Monoclonal
Clone	TRAP1-6
Preservative	0.05% Sodium Azide
Isotype	IgG1
Purity	Immunogen affinity purified
Buffer	PBS with 1 mg/ml BSA

Product Description	
Host	Mouse
Gene ID	10131
Gene Symbol	TRAP1
Species	Human, Mouse
Reactivity Notes	Mouse reported in literature (PMID: 21924824). Please note that this antibody is reactive to Mouse and derived from the same host, Mouse. Additional Mouse on Mouse blocking steps may be required for IHC and ICC experiments. Please contact Technical Support for more information.
Specificity/Sensitivity	Detects tumor necrosis factor receptor-associated protein (TRAP1) from human tissues.
Immunogen	Purified recombinant TRAP1.


Product Application Details	
Applications	Western Blot, Immunocytochemistry/ Immunofluorescence, Immunohistochemistry, Immunohistochemistry-Frozen, Immunohistochemistry-Paraffin, Immunoprecipitation
Recommended Dilutions	Western Blot 1:2000, Immunohistochemistry 1:10 - 1:500, Immunocytochemistry/ Immunofluorescence 1:250, Immunoprecipitation 1:10 - 1:500, Immunohistochemistry-Paraffin 1:10 - 1:100, Immunohistochemistry-Frozen 1:20
Application Notes	WB: Detects an approx. 75 kDa protein representing TRAP1.

Images


Western Blot: TRAP1 Antibody (TRAP1-6) [NB300-555] - Figure 2 illustrates IP of TRAP1 visualized by Coomassie Blue staining.

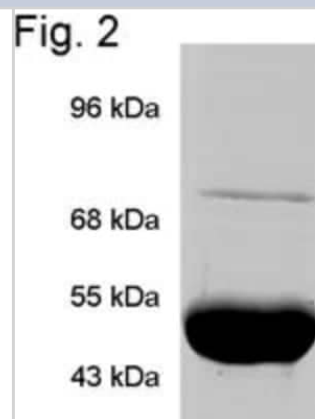

Immunocytochemistry/Immunofluorescence: TRAP1 Antibody (TRAP1-6) [NB300-555] - Analysis of TRAP1 in NCI-H460 Cells. Cells were grown on chamber slides and fixed with formaldehyde prior to staining. Cells were probed without (control) or with a TRAP1 monoclonal antibody at a dilution of 1:200 overnight at 4C, washed with PBS and incubated with a DyLight-488 conjugated secondary antibody. TRAP1 staining (green), F-Actin staining with Phalloidin (red) and nuclei with DAPI (blue) is shown.

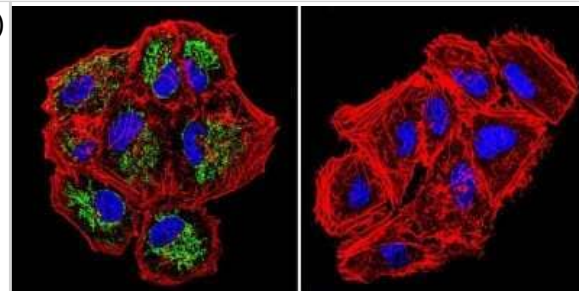


Immunohistochemistry-Paraffin: TRAP1 Antibody (TRAP1-6) [NB300-555] - Normal biopsies of deparaffinized Human liver tissue, dilution of 1:20.

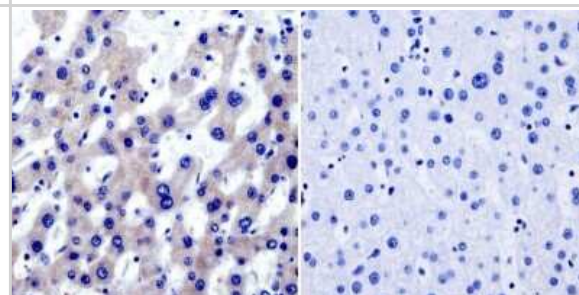


Immunocytochemistry/Immunofluorescence: TRAP1 Antibody (TRAP1-6) [NB300-555] - Immunolocalization of TRAP1 in PC-3-M cells.


Fig. 1

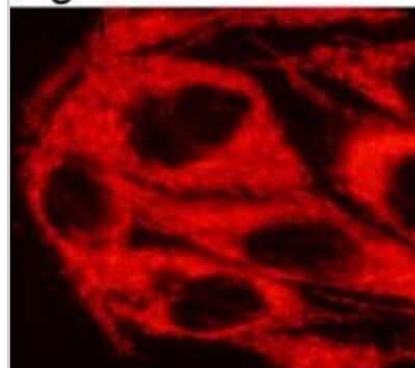

Immunohistochemistry: TRAP1 Antibody (TRAP1-6) [NB300-555] - Representative images of IHF staining of EP2 or EP4 (green) and TRAP1 (NB300-555) (red) in subchondral bone 2 weeks post-sham or ACLT surgery (left) and quantitative analysis (right). The white arrows indicate TRAP-positive osteoclasts expressing EP2 or EP4. Scale bars, 20 μ m. Image collected and cropped by CiteAb from the following publication (<https://pubmed.ncbi.nlm.nih.gov/35260562/>) licensed under a CC-BY license.

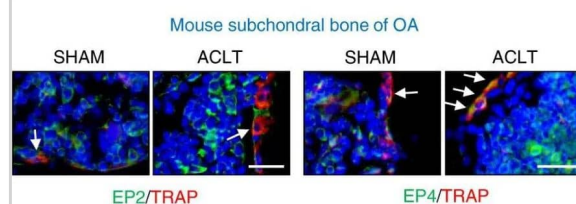

Immunocytochemistry/ Immunofluorescence: TRAP1 Antibody (TRAP1-6) [NB300-555] - PGE2 regulates migration & osteoclast differentiation through the Gαs/PI3K/AKT/MAPK signaling pathways downstream of EP4. PGE2 regulates migration & osteoclast differentiation independent of β-arrestin1 (a) & β-arrestin2 (b) but via Gαs (c). BMMs from the WT & littermate β-arrestin1- & β-arrestin2-knockout mice were used to generate osteoclasts by stimulation with 10 ng·mL⁻¹ M-CSF & 50 ng·mL⁻¹ RANKL & incubation with 100 nmol·L⁻¹ PGE2. For Gαs experiments, osteoclasts were generated from the GMMs of the Gαs^{fl/fl} controls & Gαs^{fl/fl}; LysM-cre (Gαs^{LysM}) mice. A differentiation assay (TRAP staining) (left) & the corresponding quantitative analysis (right). Scale bars, 50 μm. Error bars are the mean ± s.d. n = 3. Two-way ANOVA followed by Tukey's t tests. Scale bars, 50 μm. d cAMP production was measured via ELISAs for osteoclasts derived from the Ep4^{fl/fl} & Ep4^{LysM} mice. BMMs from the two mouse strains were isolated & stimulated with 10 ng·mL⁻¹ M-CSF & 50 ng·mL⁻¹ RANKL to differentiate into osteoclasts & incubated with 100 nmol·L⁻¹ PGE2 for 30 min prior to cAMP measurements. Error bars are the mean ± s.d. ***P < 0.001 by unpaired two-tailed Student's t test. The experiment was performed with three biological replicates. e Representative images of the indicated protein expression by western blotting for osteoclasts generated using BMMs of Ep4^{fl/fl} & Ep4^{LysM} mice. The cells were treated either with osteoclastogenic media (10 ng·mL⁻¹ M-CSF & 50 ng·mL⁻¹ RANKL) alone, with PGE2 (100 nmol·L⁻¹), or PGE2 with IBMX (1 mmol·L⁻¹) for 3 h. The experiments were performed with three biological replicates. Image collected & cropped by CiteAb from the following publication (<https://pubmed.ncbi.nlm.nih.gov/35260562>), licensed under a CC-BY license. Not internally tested by Novus Biologicals.

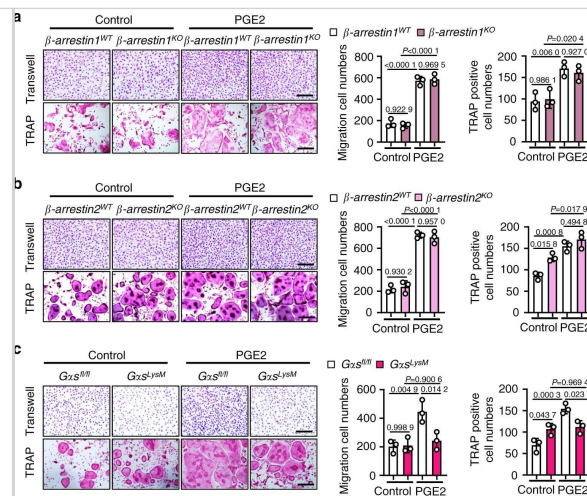

Publications

Jiang W, Jin y, Zhang S et al. PGE2 activates EP4 in subchondral bone osteoclasts to regulate osteoarthritis Bone research 2022-03-09 [PMID: 35260562] (IHC-P, Mouse)

Wang TS, Coppens I, Saorin A et al. Endolysosomal Targeting of Mitochondria Is Integral to BAX-Mediated Mitochondrial Permeabilization during Apoptosis Signaling Dev. Cell 2020-06-22 [PMID: 32504557]

Hamacher-Brady A, Brady NR et al. Bax/Bak-Dependent, Drp1-Independent Targeting of XIAP into Inner-Mitochondrial Compartments Counteracts Smac-Dependent Effector Caspase Activation J. Biol. Chem. 2015-07-01 [PMID: 26134559] (ICC/IF, Human)

Hamacher-Brady A, Choe S C, Krijnse-Locker J, Brady N R. Intramitochondrial recruitment of endolysosomes mediates Smac degradation and constitutes a novel intrinsic apoptosis antagonizing function of XIAP E3 ligase. Cell Death Differ. 2014-12-01 [PMID: 25080938]



Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112
USA
Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave
Toronto, ON M8Z 4E6
Canada
Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB, United Kingdom
Phone: (44) (0) 1235 529449
Free Phone: 0800 37 34 15
Fax: (44) (0) 1235 533420
info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com
Technical Support: nb-technical@bio-techne.com
Orders: nb-customerservice@bio-techne.com
General: novus@novusbio.com

Products Related to NB300-555

NBL1-17255	TRAP1 Overexpression Lysate
HAF007	Goat anti-Mouse IgG Secondary Antibody [HRP]
NB720-B	Rabbit anti-Mouse IgG (H+L) Secondary Antibody [Biotin]
NBP1-97005-0.5mg	Mouse IgG1 Isotype Control (MG1)

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NB300-555

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

