

# Product Datasheet

## PMEL17/SILV Antibody (P14-V)

### NBP1-79001-100ul

Unit Size: 100 ul

Store at 4C. Do not freeze.

[www.novusbio.com](http://www.novusbio.com)



[technical@novusbio.com](mailto:technical@novusbio.com)

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:  
[www.novusbio.com/NBP1-79001](http://www.novusbio.com/NBP1-79001)

Updated 7/16/2023 v.20.1

Earn rewards for product  
reviews and publications.

Submit a publication at [www.novusbio.com/publications](http://www.novusbio.com/publications)

Submit a review at [www.novusbio.com/reviews/destination/NBP1-79001](http://www.novusbio.com/reviews/destination/NBP1-79001)



**NBP1-79001-100ul**

PMEL17/SILV Antibody (P14-V)


Product Information	
Unit Size	100 ul
Concentration	Please see the vial label for concentration. If unlisted please contact technical services.
Storage	Store at 4C. Do not freeze.
Clonality	Monoclonal
Clone	P14-V
Preservative	0.05% Sodium Azide
Isotype	IgG
Purity	Immunogen affinity purified
Buffer	20mM Tris-HCl (pH 8.0) and 20mg/ml BSA

Product Description	
Host	Rabbit
Gene ID	6490
Gene Symbol	PMEL
Species	Human
Marker	Melanosome Marker
Immunogen	Peptide derived from C-terminal sequence of human PMEL17/SILV/gp100. Antibody recognizes the epitope between Ile649 - Val661.
Notes	This antibody is immunoaffinity purified with immunogenic peptide as a ligand.

Product Application Details	
Applications	Immunohistochemistry, Immunohistochemistry-Paraffin
Recommended Dilutions	Immunohistochemistry 1:100-1:200, Immunohistochemistry-Paraffin 1:100-1:200
Application Notes	Although not tested, this antibody may also be useful in Immunohistochemistry-Frozen.

**Images**

Immunohistochemistry-Paraffin: PMEL17/SILV Antibody (P14-V) [NBP1-79001] - Diffuse anti Melanosome antibody positivity in the cutaneous malignant melanoma. Formalin fixed, paraffin embedded human tissue (4 um section) stained with anti -Melanosome (DB 049) monospecific clonal antibody.

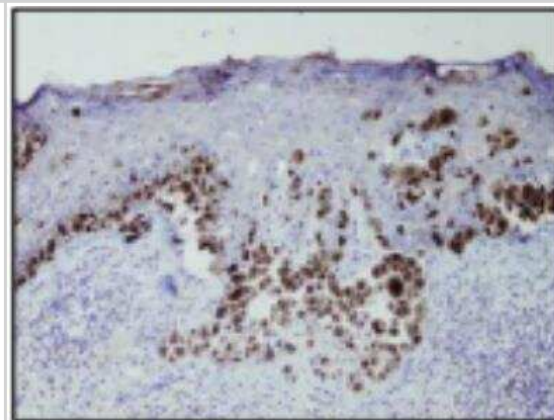

## Procedures

### Immunohistochemistry-Paraffin protocol for PMEL17/SILV Antibody (NBP1-79001)

IHC-P protocol (NBP1-79001):

1. Deparaffinize the section in 3 changes of xylene, 5 minutes each.
2. Wash the section in 96%, 80% and 70% benzyl alcohol for 5 minutes each.
3. Rinse in distilled water.
4. Block the endogenous peroxidase by incubating the tissue in 3% hydrogen peroxide (H<sub>2</sub>O<sub>2</sub>) for 10 minutes.
5. Wash in distilled water for 5 minutes.
6. For antigen retrieval: immerse the slide in the citrate buffer, pH 6.0, 0.05% Tween-20\*, and incubate in microwave (600W) for 10 minutes. (Alternatively adjust to your own protocol, keeping the required pH)
7. Remove the staining to room temperature and let the slide to cool (in citrate buffer, pH 6.0) for 15 minutes.
8. Rinse in distilled water.
9. Wash in 0.05 M Tris-HCl , pH 7.6 buffer supplemented with 0.2% of Tween-20 (buffer A) for 5 minutes.
10. Incubate the section with primary antibody diluted in buffer A\*\* at the dilution 1:100 - 200 for 1 hour in the closed wet chamber.
11. Wash twice 5 minutes with buffer A.
12. Apply the secondary antibody (the protocol depends on the supplier), and proceed to standard immunohistochemistry protocol (HRP - Peroxide - DAB).
13. Wash twice 5 minutes with buffer A.
14. Apply the chromogen (DAB), 10 minutes.
15. Wash in water for 10 minutes.
16. Stain in hematoxylin for 5 minutes.
17. Wash in water for 10 minutes.
18. Dehydrate the section in 2 changes of 96% benzyl alcohol for 5 minutes each.
19. Wash the section in 2 changes of xylene for 2 minutes each.
20. Mount the slide for observation.





### **Novus Biologicals USA**

10730 E. Briarwood Avenue  
Centennial, CO 80112  
USA  
Phone: 303.730.1950  
Toll Free: 1.888.506.6887  
Fax: 303.730.1966  
nb-customerservice@bio-techne.com

### **Bio-Techne Canada**

21 Canmotor Ave  
Toronto, ON M8Z 4E6  
Canada  
Phone: 905.827.6400  
Toll Free: 855.668.8722  
Fax: 905.827.6402  
canada.inquires@bio-techne.com

### **Bio-Techne Ltd**

19 Barton Lane  
Abingdon Science Park  
Abingdon, OX14 3NB, United Kingdom  
Phone: (44) (0) 1235 529449  
Free Phone: 0800 37 34 15  
Fax: (44) (0) 1235 533420  
info.EMEA@bio-techne.com

### **General Contact Information**

www.novusbio.com  
Technical Support: nb-technical@bio-techne.com  
Orders: nb-customerservice@bio-techne.com  
General: novus@novusbio.com

### **Limitations**

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit [www.novusbio.com/guarantee](http://www.novusbio.com/guarantee)

Earn gift cards/discounts by submitting a review: [www.novusbio.com/reviews/submit/NBP1-79001](http://www.novusbio.com/reviews/submit/NBP1-79001)

Earn gift cards/discounts by submitting a publication using this product:  
[www.novusbio.com/publications](http://www.novusbio.com/publications)

