

Product Datasheet

GPR175 Antibody NLS1522

Unit Size: 0.05 ml

Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/NLS1522

Updated 4/18/2021 v.20.1

Earn rewards for product
reviews and publications.

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/NLS1522



NLS1522**GPR175 Antibody**


Product Information	
Unit Size	0.05 ml
Concentration	1.0 mg/ml
Storage	Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.
Clonality	Polyclonal
Preservative	0.1% Sodium Azide
Isotype	IgG
Purity	Immunogen affinity purified
Buffer	PBS

Product Description	
Host	Rabbit
Gene ID	131601
Gene Symbol	TPRA1
Species	Human, Primate, Monkey
Reactivity Notes	Gorilla, Gibbon, Marmoset. Predicted cross-reactivity based on sequence identity: Elephant (94%), Seabass (94%), Salmon (94%), Pufferfish (94%), Stickleback (94%), Mouse (88%), Rat (88%), Canine (88%), Bovine (88%), Panda (88%), Rabbit (88%), Porcine (88%), Platypus (88%), Lizard (88%).
Specificity/Sensitivity	Human TPRA1 / GPR175. BLAST analysis of the peptide immunogen showed no homology with other human proteins.
Immunogen	Synthetic 16 amino acid peptide from 1st cytoplasmic domain of human TPRA1 / GPR175.

Product Application Details	
Applications	Immunohistochemistry, Immunohistochemistry-Paraffin, Immunocytochemistry
Recommended Dilutions	Immunohistochemistry, Immunohistochemistry-Paraffin 4-19 ug/ml, Immunocytochemistry 1:10 - 1:500

Images


Immunohistochemistry-Paraffin: GPR175 Antibody [NLS1522] - Anti-TPRA1 / GPR175 antibody IHC of human kidney. Immunohistochemistry of formalin-fixed, paraffin-embedded tissue after heat-induced antigen retrieval.

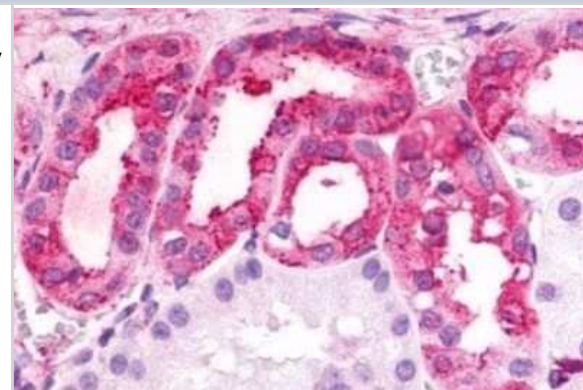


Immunohistochemistry-Paraffin: GPR175 Antibody [NLS1522] - Analysis of anti-TPRA1 / GPR175 antibody with cells expressing vector only at 8 ug/ ml.




Immunohistochemistry: GPR175 Antibody [NLS1522] - Analysis of transfected cells expressing GPR175 using this antibody at 8 ug/ml.

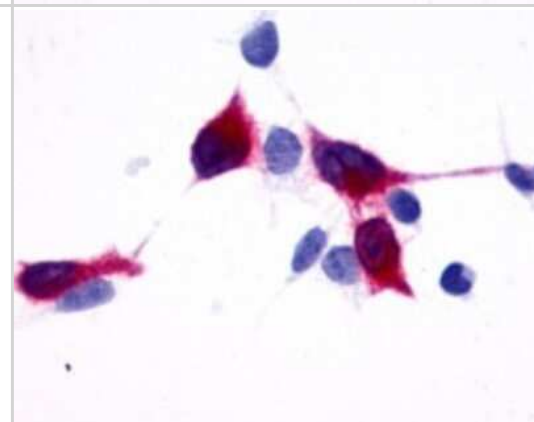


Immunohistochemistry-Paraffin: GPR175 Antibody [NLS1522] - IHC analysis of Kidney using this antibody at 8 ug/ml.

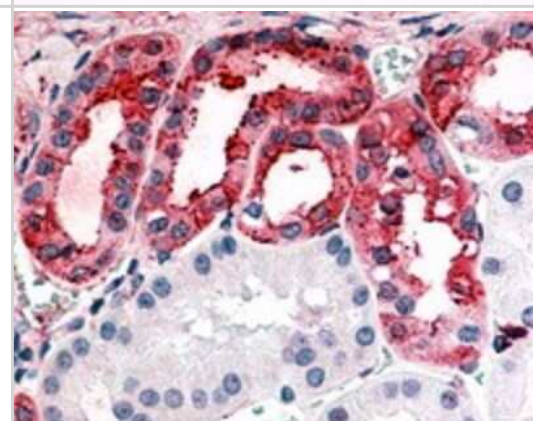


Immunohistochemistry-Paraffin: GPR175 Antibody [NLS1522] - Analysis of anti-TPRA1 / GPR175 antibody with human brain, glioblastoma.

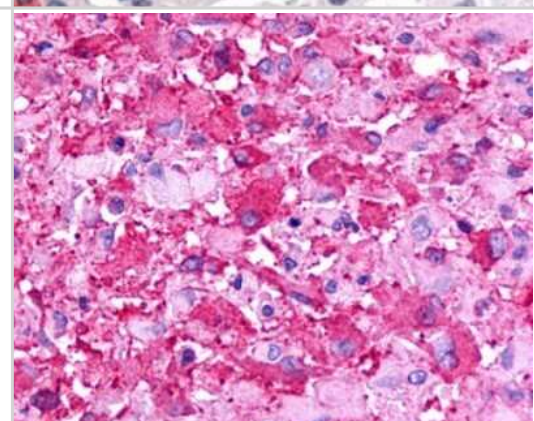


Immunohistochemistry-Paraffin: GPR175 Antibody [NLS1522] - Anti-TPRA1 / GPR175 antibody IHC of human Breast, Carcinoma. Immunohistochemistry of formalin-fixed, paraffin-embedded tissue after heat-induced antigen retrieval.

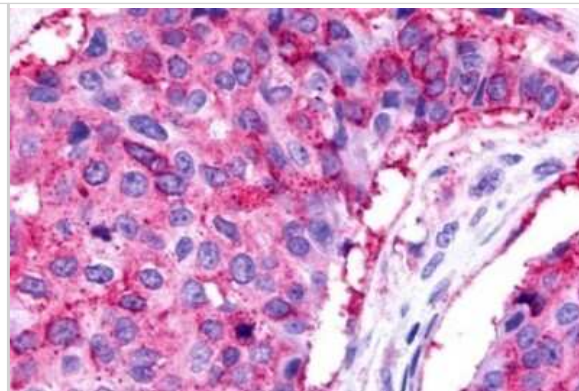


Immunohistochemistry-Paraffin: GPR175 Antibody [NLS1522] - Anti-TPRA1 / GPR175 antibody IHC of human Lung, Adenocarcinoma. Immunohistochemistry of formalin-fixed, paraffin-embedded tissue after heat-induced antigen retrieval.

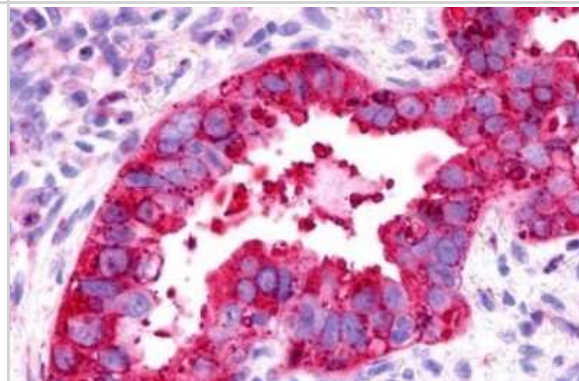

Immunohistochemistry-Paraffin: GPR175 Antibody [NLS1522] - Anti-TPRA1 / GPR175 antibody IHC of human Pancreas, Carcinoma. Immunohistochemistry of formalin-fixed, paraffin-embedded tissue after heat-induced antigen retrieval.

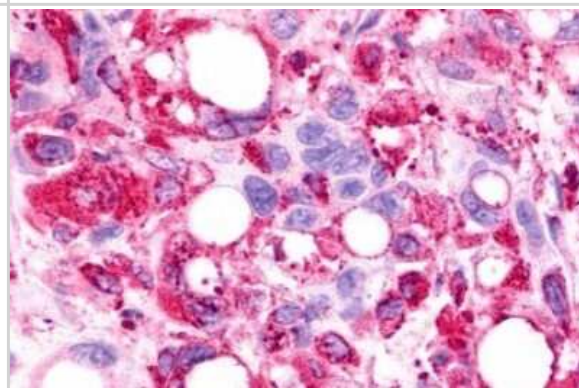

Immunohistochemistry-Paraffin: GPR175 Antibody [NLS1522] - Anti-TPRA1 / GPR175 antibody IHC of human heart. Immunohistochemistry of formalin-fixed, paraffin-embedded tissue after heat-induced antigen retrieval.

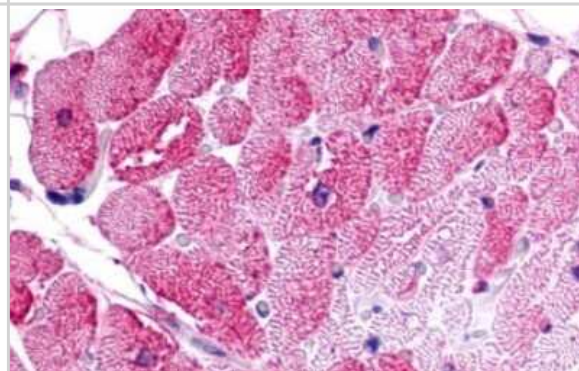

Procedures

Immunohistochemistry protocol for GPR175 Antibody (NLS1522)

Immunohistochemistry Protocol for GPR175 Antibody (NLS1522):

Immunohistochemistry

1. Prepare tissue with formalin fixation and by embedding it in paraffin wax.
2. Make 4-um sections and place on pre-cleaned and charged microscope slides.
3. Heat in a tissue-drying oven for 45 minutes at 60 degrees Celcius.
4. Deparaffinize the tissues by wash drying the slides in 3 changes of xylene approximately 5 minutes each @ RT.
5. Rehydrate the tissues by washing the slides in 3 changes of 100% alcohol approximately 3 minutes each @ RT.
6. Wash the slides in 2 changes of 95% alcohol approximately 3 minutes each @ RT.
7. Wash the slides in 1 change of 80% alcohol approximately 3 minutes @ RT.
8. Rinse the slides in gentle running distilled water approximately 5 minutes @ RT.
9. Perform antigen retrieval by steaming the slides in 0.01M sodium citrate buffer (pH 6.0) @ 99-100 degrees Celcius for 20 minutes.
10. Remove the slides from the heat and let stand in buffer @ RT for 20 minutes.
11. Rinse the slides in 1X TBS-T for 1 minute @ RT.

****Do not allow the tissues to dry at any time during the staining procedure****

12. Begin the immunostaining by applying a universal protein block approximately 20 minutes @ RT.
13. Drain protein block from the slides and apply the diluted primary antibody approximately 45 minutes @ RT.
14. Rinse the slide in 1X TBS-T approximately 1 minute @ RT.
15. Apply a biotinylated anti-rabbit IgG (H+L) secondary approximately 30 minutes @ RT.
16. Rinse the slide in 1X TBS-T approximately 1 minute at RT.
17. Apply an alkaline phosphatase streptavidin approximately 30 minutes at RT.
18. Rinse the slide in 1X TBS-T approximately 1 minute at RT.
19. Apply an alkaline phosphatase chromagen substrate approximately 30 minutes at RT.
20. Rinse the slide in distilled water approximately 1 minute @ RT.

****This method should only be used if the chromagen substrate is alcohol insoluble (ie: Vector Red, DAB)****

21. Dehydrate the tissue by washing the slides in 2 changes of 80% alcohol approximately 1 minute each @ RT.
22. Wash the slides in 2 changes of 95% alcohol approximately 1 minute each @ RT.
23. Wash the slides in 3 changes of 100% alcohol approximately 1 minute each @ RT.
24. Wash the slides in 3 changes of xylene approximately 1 minute each @ RT.
25. Apply cover slip.





Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112
USA
Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave
Toronto, ON M8Z 4E6
Canada
Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB, United Kingdom
Phone: (44) (0) 1235 529449
Free Phone: 0800 37 34 15
Fax: (44) (0) 1235 533420
info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com
Technical Support: nb-technical@bio-
techne.com
Orders: nb-customerservice@bio-techne.com
General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NLS1522

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

