Product Datasheet

Leukotriene B4 Receptor 2 Antibody NLS2099

Unit Size: 0.05 ml

Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at: www.novusbio.com/NLS2099

Updated 7/26/2021 v.20.1

Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications
Submit a review at www.novusbio.com/reviews/destination/NLS2099




NLS2099

Leukotriene B4 Receptor 2 Antibody

<u>'</u>	•
Product Information	
Unit Size	0.05 ml
Concentration	1.0 mg/ml
Storage	Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.
Clonality	Polyclonal
Preservative	0.1% Sodium Azide
Isotype	IgG
Purity	Immunogen affinity purified
Buffer	PBS
Product Description	
Host	Rabbit
Gene ID	56413
Gene Symbol	LTB4R2
Species	Human, Monkey, Primate
Reactivity Notes	Predicted cross-reactivity based on sequence identity: Gorilla (100%), Marmoset (100%), Mouse (88%), Rat (88%), Bat (88%).
Specificity/Sensitivity	Human LTB4R2. BLAST analysis of the peptide immunogen showed no homology with other human proteins.
Immunogen	Synthetic 16 amino acid peptide from 3rd cytoplasmic domain of human Leukotriene B4 Receptor 2.
Product Application Details	
Applications	Flow Cytometry, Immunohistochemistry, Immunohistochemistry-Frozen, Immunohistochemistry-Paraffin
Recommended Dilutions	Flow Cytometry 1:10 - 1:1000, Immunohistochemistry, Immunohistochemistry-Paraffin 26 ug/ml, Immunohistochemistry-Frozen 1:10 - 1:500


Images

Immunohistochemistry-Paraffin: Leukotriene B4 Receptor 2 Antibody [NLS2099] - Lung, Adenocarcinoma

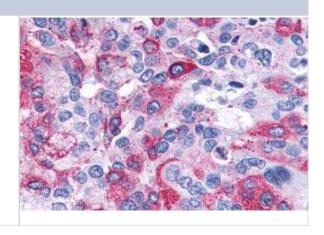





Immunohistochemistry-Paraffin: Leukotriene B4 Receptor 2 Antibody [NLS2099] - Analysis of anti-LTB4R2 / BLT2 antibody with human colon, carcinoma and leukocytes.

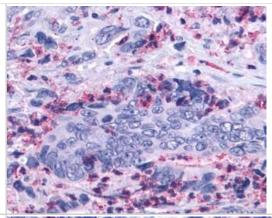


Immunohistochemistry-Paraffin: Leukotriene B4 Receptor 2 Antibody [NLS2099] - Human Lymph Node, Hodgkins Lymphoma. Immunohistochemistry of formalin-fixed, paraffin-embedded tissue after heat-induced antigen retrieval.

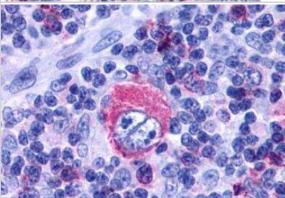


Immunohistochemistry-Paraffin: Leukotriene B4 Receptor 2 Antibody [NLS2099] - Anti-LTB4R2 antibody IHC of human tonsil. Immunohistochemistry of formalin-fixed, paraffin-embedded tissue after heat-induced antigen retrieval.

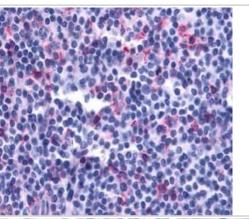

Immunohistochemistry-Paraffin: Leukotriene B4 Receptor 2 Antibody [NLS2099] - Breast, carcinoma

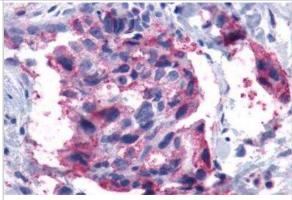

Immunohistochemistry-Paraffin: Leukotriene B4 Receptor 2 Antibody [NLS2099] - Colon, carcinoma

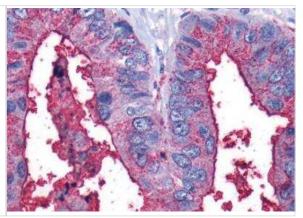

Procedures

Immunohistochemistry protocol for Leukotriene B4 Receptor 2 Antibody (NLS2099)

Immunohistochemistry Protocol for Leukotriene B4 Receptor 2 Antibody (NLS2099): Immunohistochemistry

- 1. Prepare tissue with formalin fixation and by embedding it in paraffin wax.
- 2. Make 4-um sections and place on pre-cleaned and charged microscope slides.
- 3. Heat in a tissue-drying oven for 45 minutes at 60 degrees Celcius.
- 4. Deparaffinize the tissues by wash drying the slides in 3 changes of xylene approximately 5 minutes each @ RT.
- 5. Rehydrate the tissues by washing the slides in 3 changes of 100% alcohol approximately 3 minutes each @ RT.
- Wash the slides in 2 changes of 95% alcohol approximately 3 minutes each @ RT.
- 7. Wash the slides in 1 change of 80% alcohol approximately 3 minutes @ RT.
- 8. Rinse the slides in gentle running distilled water approximately 5 minutes @ RT.
- 9. Perform antigen retrieval by steaming the slides in 0.01M sodium citrate buffer (pH 6.0) @ 99-100 degrees Celcius for 20 minutes.
- 10. Remove the slides from the heat and let stand in buffer @ RT for 20 minutes.
- 11. Rinse the slides in 1X TBS-T for 1 minute @ RT.
- **Do not allow the tissues to dry at any time during the staining procedure**
- 12. Begin the immunostaining by applying a universal protein block approximately 20 minutes @ RT.
- 13. Drain protein block from the slides and apply the diluted primary antibody approximately 45 minutes @ RT.
- 14. Rinse the slide in 1X TBS-T approximately 1 minute @ RT.
- 15. Apply a biotinylated anti-rabbit IgG (H+L) secondary approximately 30 minutes @ RT.
- 16. Rinse the slide in 1X TBS-T approximately 1 minute at RT.
- 17. Apply an alkaline phosphatase steptavidin approximately 30 minutes at RT.
- 18. Rinse the slide in 1X TBS-T approximately 1 minute at RT.
- 19. Apply an alkaline phosphatase chromagen substrate approximately 30 minutes at RT.
- 20. Rinse the slide in distilled water approximately 1 minute @ RT.
- **This method should only be used if the chromagen substrate is alcohol insoluble (ie: Vector Red, DAB)**
- 21. Dehydrate the tissue by washing the slides in 2 changes of 80% alcohol approximately 1 minute each @ RT.
- 22. Wash the slides in 2 changes of 95% alcohol approximately 1 minute each @ RT.
- 23. Wash the slides in 3 changes of 100% alcohol approximately 1 minute each @ RT.
- 24. Wash the slides in 3 changes of xyleneapproximately 1 minute each @ RT.
- 25. Apply cover slip.





Novus Biologicals USA

10730 E. Briarwood Avenue Centennial, CO 80112 USA

Phone: 303.730.1950 Toll Free: 1.888.506.6887

Fax: 303.730.1966

nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave Toronto, ON M8Z 4E6 Canada

Phone: 905.827.6400 Toll Free: 855.668.8722 Fax: 905.827.6402

canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane Abingdon Science Park Abingdon, OX14 3NB, United Kingdom Phone: (44) (0) 1235 529449

Free Phone: 0800 37 34 15 Fax: (44) (0) 1235 533420 info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com Technical Support: nb-technical@biotechne.com

Orders: nb-customerservice@bio-techne.com

General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NLS2099

Earn gift cards/discounts by submitting a publication using this product: www.novusbio.com/publications

