Product Datasheet

SF-1/NR5A1/Steroidogenic Factor 1 Antibody - BSA Free NLS5534

Unit Size: 0.05 ml

Keep as concentrated solution. Aliquot and store at -20C or below. Avoid multiple freeze-thaw cycles.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at: www.novusbio.com/NLS5534

Updated 10/23/2024 v.20.1

Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications
Submit a review at www.novusbio.com/reviews/destination/NLS5534



NLS5534


SF-1/NR5A1/Steroidogenic Factor 1 Antibody - BSA Free		
Product Information		
Product Description		
on Y		

Product Application Details		
Applications	Immunohistochemistry, Immunohistochemistry-Paraffin	
Recommended Dilutions	Immunohistochemistry, Immunohistochemistry-Paraffin 1:100	




Images

Immunohistochemistry-Paraffin: SF-1/NR5A1/Steroidogenic Factor 1 Antibody [NLS5534] - Analysis of anti-NR5A1 antibody with human ovary.

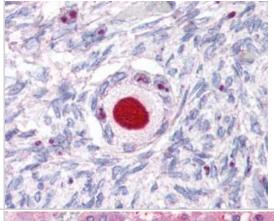

Immunohistochemistry-Paraffin: SF-1/NR5A1/Steroidogenic Factor 1 Antibody - BSA Free [NLS5534] - Immunohistochemistry of formalin-fixed, paraffin-embedded Human Adrenal tissue

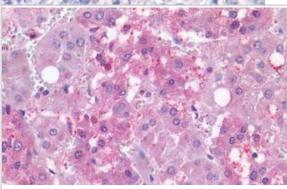

Procedures

Immunohistochemistry Protocol for Steroidogenic Factor 1 Antibody (NLS5534)

Immunohistochemistry Protocol for Steroidogenic Factor 1 Antibody (NLS5534): Immunohistochemistry

- 1. Prepare tissue with formalin fixation and by embedding it in paraffin wax.
- 2. Make 4-um sections and place on pre-cleaned and charged microscope slides.
- 3. Heat in a tissue-drying oven for 45 minutes at 60 degrees Celcius.
- 4. Deparaffinize the tissues by wash drying the slides in 3 changes of xylene approximately 5 minutes each @ RT.
- 5. Rehydrate the tissues by washing the slides in 3 changes of 100% alcohol approximately 3 minutes each @ RT.
- 6. Wash the slides in 2 changes of 95% alcohol approximately 3 minutes each @ RT.
- 7. Wash the slides in 1 change of 80% alcohol approximately 3 minutes @ RT.
- 8. Rinse the slides in gentle running distilled water approximately 5 minutes @ RT.
- 9. Perform antigen retrieval by steaming the slides in 0.01M sodium citrate buffer (pH 6.0) @ 99-100 degrees Celcius for 20 minutes.
- 10. Remove the slides from the heat and let stand in buffer @ RT for 20 minutes.
- 11. Rinse the slides in 1X TBS-T for 1 minute @ RT.
- **Do not allow the tissues to dry at any time during the staining procedure**
- 12. Begin the immunostaining by applying a universal protein block approximately 20 minutes @ RT.
- 13. Drain protein block from the slides and apply the diluted primary antibody approximately 45 minutes @ RT.
- 14. Rinse the slide in 1X TBS-T approximately 1 minute @ RT.
- 15. Apply a biotinylated anti-rabbit IgG (H+L) secondary approximately 30 minutes @ RT.
- 16. Rinse the slide in 1X TBS-T approximately 1 minute at RT.
- 17. Apply an alkaline phosphatase steptavidin approximately 30 minutes at RT.
- 18. Rinse the slide in 1X TBS-T approximately 1 minute at RT.
- 19. Apply an alkaline phosphatase chromagen substrate approximately 30 minutes at RT.
- 20. Rinse the slide in distilled water approximately 1 minute @ RT.
- **This method should only be used if the chromagen substrate is alcohol insoluble (ie: Vector Red, DAB)**
- 21. Dehydrate the tissue by washing the slides in 2 changes of 80% alcohol approximately 1 minute each @ RT.
- 22. Wash the slides in 2 changes of 95% alcohol approximately 1 minute each @ RT.
- 23. Wash the slides in 3 changes of 100% alcohol approximately 1 minute each @ RT.
- 24. Wash the slides in 3 changes of xyleneapproximately 1 minute each @ RT.
- 25. Apply cover slip.





Novus Biologicals USA

10730 E. Briarwood Avenue Centennial, CO 80112

USA

Phone: 303.730.1950 Toll Free: 1.888.506.6887

Fax: 303.730.1966

nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave Toronto, ON M8Z 4E6

Canada

Phone: 905.827.6400 Toll Free: 855.668.8722 Fax: 905.827.6402

canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane Abingdon Science Park Abingdon, OX14 3NB, United Kingdom Phone: (44) (0) 1235 529449

Free Phone: 0800 37 34 15 Fax: (44) (0) 1235 533420 info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com

Technical Support: nb-technical@bio-

techne.com

Orders: nb-customerservice@bio-techne.com

General: novus@novusbio.com

Products Related to NLS5534

HAF008 Goat anti-Rabbit IgG Secondary Antibody [HRP]

NB7160 Goat anti-Rabbit IgG (H+L) Secondary Antibody [HRP]

NBP2-24891 Rabbit IgG Isotype Control

H00002516-P01-10ug Recombinant Human SF-1/NR5A1/Steroidogenic Factor 1 GST (N-

Term) Protein

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NLS5534

Earn gift cards/discounts by submitting a publication using this product: www.novusbio.com/publications

