Product Datasheet

Lightning-Link (R) Rapid AMCA Antibody Labeling Kit 313-0015

Unit Size: 1 mg Store at -20C.

www.novusbio.com



technical@novusbio.com

Publications: 1

Protocols, Publications, Related Products, Reviews, Research Tools and Images at: www.novusbio.com/313-0015

Updated 5/23/2024 v.20.1

Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications
Submit a review at www.novusbio.com/reviews/destination/313-0015



313-0015

Lightning-Link (R) Rapid A	AMCA Antibody Labeling Kit
Product Information	
Unit Size	1 mg
Concentration	Concentration is not relevant for this product. Please see the protocols for proper use of this product.
Storage	Store at -20C.
Conjugate	AMCA
Product Description	
Description	Lightning-Link Rapid is an innovative technology that enables direct labeling of proteins, peptides or other biomolecules for use in R&D applications, drug discovery and the development of diagnostic kits (See protocol for further information). The easy-to-use, one step procedure allows researchers to covalently label biomolecules with only 30 seconds hands-on time; furthermore conjugates are ready to use in less than twenty minutes. The researcher simply pipettes the biomolecule into a vial of lyophilized mixture containing the label of interest and incubates (for more details please watch the video below). FeaturesBenefitsQuick and easy to useSave time, no special knowledge requiredNo separation steps100% recovery - no antibody/protein lossCan be used in a wide range of applicationsFlexibleFreeze driedShips at ambient temperature, long shelf-lifeFully scalable (10 ug to 1 g or more)Easy transfer from R&D to manufacturingStringently QC testedConsistent high quality, excellent batch-to-batch reproducibilityLarge number of labels available Experimental flexibilityReliable: nearly 300 referencesSuccessfully used in many fields of research Unlike standard AMCA labeling procedures, where a large molar excess of AMCA is employed, the Lightning-Link AMCA kit combines a low molar ratio of AMCA:antibody and high labeling efficiency to eliminate the need for purification steps. AMCA (aminomethylcoumarin acetate) is a blue fluorescent dye commonly used in fluorescence microscopy, arrays and in-situ hybridisation. It has an excitation wavelength of 352nm and a maximal emission of 452nm. Learn more about Lightning-Link™ Conjugation Kits by reading FAQs For more information please check out these useful links! Antibody Labeling Guide
	Antibody Conjugation Illustrated Assay
Kit Components	1 or 3 glass vial(s) of Lightning-Link Rapid mix, 1 vial of LL-Rapid Modifier reagent, 1 vial of LL-Rapid Quencher reagent



Notes

This product is manufactured by Abcam and distributed by Novus Biologicals.

This product is for research use only and is not approved for use in humans or in clinical diagnosis. This product is guaranteed for 1 year from date of receipt and this statement overrides any mentioned guarantee period on the limitations section of this products datasheet. Please contact technical@novusbio.com with questions.


Product Application Details

Application Notes

By circumventing the desalting or dialysis steps that commonly interrupt traditional antibody conjugation procedures, LightningLink technology can be used to label both small (e.g. 10 ug) and large quantities of primary antibodies with ease. Batch-to-batch variation upon scale up is minimal as the process is so simple, and recoveries are always 100%. This kit can be used to label up to 2mg of antibody, and is supplied in one vial.

Images

Lightning-Link Rapid AMCA Antibody Labeling Kit [313-0015] - Lightning-Link Rapid AMCA Antibody Labeling Kit [313-0010] - AMCA (aminomethylcoumarin acetate) is a blue fluorescent dye commonly used in fluorescence microscopy, arrays and in-situ hybridisation. It has a strong absorption at 353nm and high fluorescence at 442nm (extinction coefficient 1.9 x104 cm-1M-1).

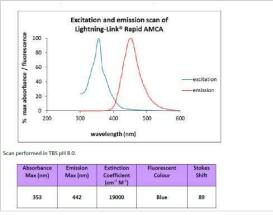

Publications

Liu Y, Wei J, Lu J et al. Micropatterned coculture of hepatocytes on electrospun fibers as a potential in vitro model for predictive drug metabolism Mater Sci Eng C Mater Biol Appl. 2016-06-01 [PMID: 27040241] (ICC/IF)





Novus Biologicals USA

10730 E. Briarwood Avenue Centennial, CO 80112 USA

Phone: 303.730.1950 Toll Free: 1.888.506.6887

Fax: 303.730.1966

nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave Toronto, ON M8Z 4E6 Canada

Phone: 905.827.6400 Toll Free: 855.668.8722 Fax: 905.827.6402

canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane Abingdon Science Park Abingdon, OX14 3NB, United Kingdom Phone: (44) (0) 1235 529449

Free Phone: 0800 37 34 15 Fax: (44) (0) 1235 533420 info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com Technical Support: nb-technical@biotechne.com

Orders: nb-customerservice@bio-techne.com

General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Kits are guaranteed for 6 months from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/313-0015

Earn gift cards/discounts by submitting a publication using this product: www.novusbio.com/publications

