

Product Datasheet

Lightning-Link (R) Rapid AMCA Antibody Labeling Kit 313-0030

Unit Size: 3 x 10ug Reaction

Store at -20C.

www.novusbio.com



technical@novusbio.com

Publications: 1

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/313-0030

Updated 5/23/2024 v.20.1

Earn rewards for product
reviews and publications.

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/313-0030



313-0030**Lightning-Link (R) Rapid AMCA Antibody Labeling Kit**

Product Information	
Unit Size	3 x 10ug Reaction
Concentration	Concentration is not relevant for this product. Please see the protocols for proper use of this product.
Storage	Store at -20C.
Conjugate	AMCA
Product Description	
Description	<p>Lightning-Link Rapid is an innovative technology that enables direct labeling of proteins, peptides or other biomolecules for use in R&D applications, drug discovery and the development of diagnostic kits (See protocol for further information).</p> <p>The easy-to-use, one step procedure allows researchers to covalently label biomolecules with only 30 seconds hands-on time; furthermore conjugates are ready to use in less than twenty minutes.</p> <p>The researcher simply pipettes the biomolecule into a vial of lyophilized mixture containing the label of interest and incubates (for more details please watch the video below).</p> <p>Features Quick and easy to use Benefits Save time, no special knowledge required No separation steps 100% recovery - no antibody/protein loss Can be used in a wide range of applications Flexible Freeze dried Ships at ambient temperature, long shelf-life Fully scalable (10 ug to 1 g or more) Easy transfer from R&D to manufacturing Stringently QC tested Consistent high quality, excellent batch-to-batch reproducibility Large number of labels available Experimental flexibility Reliable: nearly 300 references Successfully used in many fields of research</p> <p>Unlike standard AMCA labeling procedures, where a large molar excess of AMCA is employed, the Lightning-Link AMCA kit combines a low molar ratio of AMCA:antibody and high labeling efficiency to eliminate the need for purification steps.</p> <p>AMCA (aminomethylcoumarin acetate) is a blue fluorescent dye commonly used in fluorescence microscopy, arrays and in-situ hybridisation. It has an excitation wavelength of 352nm and a maximal emission of 452nm.</p> <p>Learn more about Lightning-Link™ Conjugation Kits by reading FAQs</p> <p>For more information please check out these useful links! Antibody Labeling Guide Antibody Conjugation Illustrated Assay</p>
Kit Components	1 or 3 glass vial(s) of Lightning-Link Rapid mix, 1 vial of LL-Rapid Modifier reagent, 1 vial of LL-Rapid Quencher reagent




Notes	<p>This product is manufactured by Abcam and distributed by Novus Biologicals.</p> <p>This product is for research use only and is not approved for use in humans or in clinical diagnosis. This product is guaranteed for 1 year from date of receipt and this statement overrides any mentioned guarantee period on the limitations section of this products datasheet. Please contact technical@novusbio.com with questions.</p>
--------------	---

Product Application Details

Application Notes	<p>By circumventing the desalting or dialysis steps that commonly interrupt traditional antibody conjugation procedures, LightningLink technology can be used to label both small (e.g. 10 ug) and large quantities of primary antibodies with ease. Batch-to-batch variation upon scale up is minimal as the process is so simple, and recoveries are always 100%. This kit is supplied with 3 vials.</p>
--------------------------	--

Images

Lightning-Link Rapid AMCA Antibody Labeling Kit [313-0030] - Lightning-Link Rapid AMCA Antibody Labeling Kit [313-0010] - AMCA (aminomethylcoumarin acetate) is a blue fluorescent dye commonly used in fluorescence microscopy, arrays and in-situ hybridisation. It has a strong absorption at 353nm and high fluorescence at 442nm (extinction coefficient $1.9 \times 10^4 \text{ cm}^{-1}\text{M}^{-1}$).

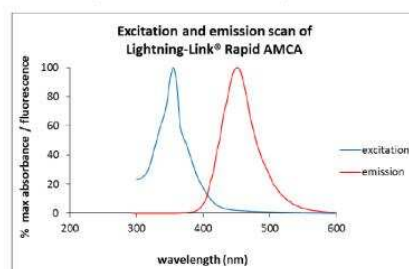

Scan performed in TBS pH 8.0.

Absorbance Max (nm)	Emission Max (nm)	Extinction Coefficient ($\text{cm}^{-1}\text{M}^{-1}$)	Fluorescent Colour	Stokes Shift
353	442	19000	Blue	89

Publications

Liu Y, Wei J, Lu J et al. Micropatterned coculture of hepatocytes on electrospun fibers as a potential in vitro model for predictive drug metabolism Mater Sci Eng C Mater Biol Appl. 2016-06-01 [PMID: 27040241] (ICC/IF)



Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112
USA
Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave
Toronto, ON M8Z 4E6
Canada
Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB, United Kingdom
Phone: (44) (0) 1235 529449
Free Phone: 0800 37 34 15
Fax: (44) (0) 1235 533420
info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com
Technical Support: nb-technical@bio-techne.com
Orders: nb-customerservice@bio-techne.com
General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Kits are guaranteed for 6 months from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/313-0030

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

