

Product Datasheet

Lightning-Link (R) Rapid Atto390 Antibody Labeling Kit 349-0010

Unit Size: 3 x 100ug Reaction

Store at -20C.

www.novusbio.com



technical@novusbio.com

Publications: 2

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/349-0010

Updated 10/23/2024 v.20.1

Earn rewards for product
reviews and publications.

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/349-0010



349-0010**Lightning-Link (R) Rapid Atto390 Antibody Labeling Kit**

Product Information	
Unit Size	3 x 100ug Reaction
Concentration	Concentration is not relevant for this product. Please see the protocols for proper use of this product.
Storage	Store at -20C.
Conjugate	Atto390
Product Description	
Description	<p>Lightning-Link Rapid is an innovative technology that enables direct labeling of proteins, peptides or other biomolecules for use in R&D applications, drug discovery and the development of diagnostic kits (See protocol for further information).</p> <p>The easy-to-use, one step procedure allows researchers to covalently label biomolecules with only 30 seconds hands-on time; furthermore conjugates are ready to use in less than twenty minutes.</p> <p>The researcher simply pipettes the biomolecule into a vial of lyophilized mixture containing the label of interest and incubates (for more details please watch the video below).</p> <p>FeaturesQuick and easy to use BenefitsSave time, no special knowledge required No separation steps 100% recovery - no antibody/protein loss Can be used in a wide range of applications Flexible Freeze dried Ships at ambient temperature, long shelf-life Fully scalable (10 ug to 1 g or more) Easy transfer from R&D to manufacturing Stringently QC tested Consistent high quality, excellent batch-to-batch reproducibility Large number of labels available Experimental flexibility Reliable: nearly 300 references Successfully used in many fields of research</p> <p>Atto390 is one of a new generation of fluorescent labels with a coumarin structure. It has a strong absorption at 388nm, high fluorescence at 468nm (extinction coefficient $2.4 \times 10^4 \text{ cm}^{-1}\text{M}^{-1}$) and high quantum yield.</p> <p>Learn more about Lightning-Link™ Conjugation Kits by reading FAQs</p> <p>For more information please check out these useful links! Antibody Labeling Guide Antibody Conjugation Illustrated Assay</p>
Kit Components	1 or 3 glass vial(s) of Lightning-Link Rapid mix, 1 vial of LL-Rapid Modifier reagent, 1 vial of LL-Rapid Quencher reagent
Notes	<p>This product is manufactured by Abcam and distributed by Novus Biologicals.</p> <p>This product is for research use only and is not approved for use in humans or in clinical diagnosis. This product is guaranteed for 1 year from date of receipt and this statement overrides any mentioned guarantee period on the limitations section of this products datasheet. Please contact technical@novusbio.com with questions.</p>
Product Application Details	




Application Notes

By circumventing the desalting or dialysis steps that commonly interrupt traditional antibody conjugation procedures, LightningLink technology can be used to label both small (e.g. 10 ug) and large quantities of primary antibodies with ease. Batch-to-batch variation upon scale up is minimal as the process is so simple, and recoveries are always 100%.

Images

Lightning-Link Rapid Atto390 Antibody Labeling Kit [349-0010]

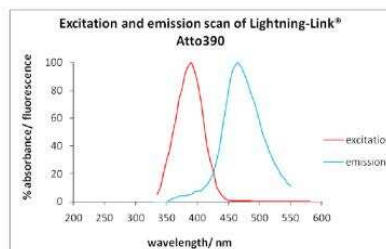

Scan performed in TBS pH 8.0.

Absorbance Max (nm)	Emission Max (nm)	Extinction Coefficient (cm ² M ⁻¹)	Fluorescent Colour	Stokes Shift
390	479	24000	Blue	89

Publications

Henning AL, McFarlin BK. Consumption of a high-fat, high-calorie meal is associated with an increase in intracellular co-localization of PPAR-gamma mRNA and protein in monocytes Methods. 2016-07-12 [PMID: 27418556] (FLOW)

Taichman RS, Patel LR, Bedenis R et al. GAS6 Receptor Status Is Associated with Dormancy and Bone Metastatic Tumor Formation. PLoS One 2013-01-01 [PMID: 23637920] (FLOW)



Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112
USA
Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave
Toronto, ON M8Z 4E6
Canada
Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB, United Kingdom
Phone: (44) (0) 1235 529449
Free Phone: 0800 37 34 15
Fax: (44) (0) 1235 533420
info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com
Technical Support: nb-technical@bio-techne.com
Orders: nb-customerservice@bio-techne.com
General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Kits are guaranteed for 6 months from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/349-0010

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

