

# Product Datasheet

## Lightning-Link (R) Rapid FluoProbes647H Antibody Labeling Kit 362-0010

Unit Size: 3 x 10ug Reaction

Store at -20C.

[www.novusbio.com](http://www.novusbio.com)



[technical@novusbio.com](mailto:technical@novusbio.com)

### Publications: 2

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:  
[www.novusbio.com/362-0010](http://www.novusbio.com/362-0010)

Updated 10/23/2024 v.20.1

Earn rewards for product  
reviews and publications.

Submit a publication at [www.novusbio.com/publications](http://www.novusbio.com/publications)

Submit a review at [www.novusbio.com/reviews/destination/362-0010](http://www.novusbio.com/reviews/destination/362-0010)



**362-0010****Lightning-Link (R) Rapid FluoProbes647H Antibody Labeling Kit**

<b>Product Information</b>	
<b>Unit Size</b>	3 x 10ug Reaction
<b>Concentration</b>	Concentration is not relevant for this product. Please see the protocols for proper use of this product.
<b>Storage</b>	Store at -20C.
<b>Conjugate</b>	FluoProbes 647H
<b>Product Description</b>	
<b>Description</b>	<p>Lightning-Link Rapid is an innovative technology that enables direct labeling of proteins, peptides or other biomolecules for use in R&amp;D applications, drug discovery and the development of diagnostic kits (See protocol for further information).</p> <p>The easy-to-use, one step procedure allows researchers to covalently label biomolecules with only 30 seconds hands-on time; furthermore conjugates are ready to use in less than twenty minutes.</p> <p>The researcher simply pipettes the biomolecule into a vial of lyophilized mixture containing the label of interest and incubates (for more details please watch the video below).</p> <p><b>Features</b> Quick and easy to use  <b>Benefits</b> Save time, no special knowledge required            No separation steps            100% recovery - no antibody/protein loss            Can be used in a wide range of applications            Flexible            Freeze dried            Ships at ambient temperature, long shelf-life            Fully scalable (10 ug to 1 g or more)            Easy transfer from R&amp;D to manufacturing            Stringently QC tested            Consistent high quality, excellent batch-to-batch reproducibility            Large number of labels available            Experimental flexibility            Reliable: nearly 300 references            Successfully used in many fields of research</p> <p>FluoProbes647H is one of a new generation of fluorescent labels designed to label biomolecules. It has a strong absorption at 650nm, high fluorescence at 684nm (extinction coefficient <math>2.5 \times 10^5 \text{ cm}^{-1}\text{M}^{-1}</math>) and high quantum yield.</p> <p>Learn more about Lightning-Link™ Conjugation Kits by reading <a href="#">FAQs</a></p> <p>For more information please check out these useful links!  <a href="#">Antibody Labeling Guide</a>  <a href="#">Antibody Conjugation Illustrated Assay</a></p>
<b>Kit Components</b>	1 or 3 glass vial(s) of Lightning-Link Rapid mix, 1 vial of LL-Rapid Modifier reagent, 1 vial of LL-Rapid Quencher reagent
<b>Notes</b>	<p>This product is manufactured by Abcam and distributed by Novus Biologicals.</p> <p>This product is for research use only and is not approved for use in humans or in clinical diagnosis. This product is guaranteed for 1 year from date of receipt and this statement overrides any mentioned guarantee period on the limitations section of this products datasheet. Please contact <a href="mailto:technical@novusbio.com">technical@novusbio.com</a> with questions.</p>


**Product Application Details**

**Application Notes**


By circumventing the desalting or dialysis steps that commonly interrupt traditional antibody conjugation procedures, LightningLink technology can be used to label both small (e.g. 10 ug) and large quantities of primary antibodies with ease. Batch-to-batch variation upon scale up is minimal as the process is so simple, and recoveries are always 100%. This kit is supplied with 3 vials, each suitable for labeling up to 200 ug of antibody.

**Images**

Flow Cytometry: Lightning-Link Rapid FluoProbes 647H Antibody Labeling Kit [362-0010] - Mouse anti-human CD3 was conjugated with FluoProbes® 647H using Lightning-Link® Rapid kit. The conjugated antibody was then used to stain human peripheral blood lymphocytes, followed by analysis with flow cytometry. (Blue line - negative control; red line - positive staining).

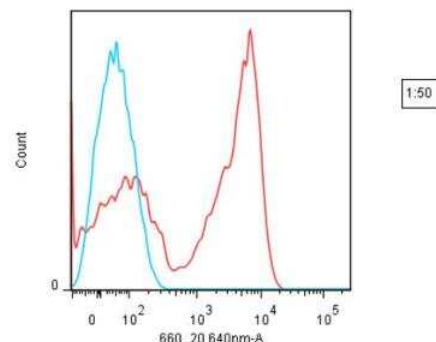

Lightning-Link Rapid FluoProbes 647H Antibody Labeling Kit [362-0010]

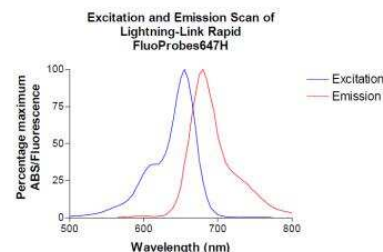

Scans performed in TBS pH 8

Absorbance Max (nm)	Emission Max (nm)	Extinction Coefficient (cm <sup>-1</sup> M <sup>-1</sup> )	Fluorescent Colour	Stokes Shift
653	674	250000	Red	21

**Publications**

Aravantinou M, Frank I, Arrode-Bruses G et al. A model of genital herpes simplex virus Type 1 infection in Rhesus Macaques. J Med Primatol. 2017-01-01 [PMID: 28748667]

Kakhlon O, Glickstein H, Feinstein N et al. Polyglucosan neurotoxicity caused by glycogen branching enzyme deficiency can be reversed by inhibition of glycogen synthase. J Neurochem 2013-04-18 [PMID: 23607684]



### **Novus Biologicals USA**

10730 E. Briarwood Avenue  
Centennial, CO 80112  
USA  
Phone: 303.730.1950  
Toll Free: 1.888.506.6887  
Fax: 303.730.1966  
nb-customerservice@bio-techne.com

### **Bio-Techne Canada**

21 Canmotor Ave  
Toronto, ON M8Z 4E6  
Canada  
Phone: 905.827.6400  
Toll Free: 855.668.8722  
Fax: 905.827.6402  
canada.inquires@bio-techne.com

### **Bio-Techne Ltd**

19 Barton Lane  
Abingdon Science Park  
Abingdon, OX14 3NB, United Kingdom  
Phone: (44) (0) 1235 529449  
Free Phone: 0800 37 34 15  
Fax: (44) (0) 1235 533420  
info.EMEA@bio-techne.com

### **General Contact Information**

www.novusbio.com  
Technical Support: nb-technical@bio-  
techne.com  
Orders: nb-customerservice@bio-techne.com  
General: novus@novusbio.com

### **Limitations**

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Kits are guaranteed for 6 months from date of receipt.

For more information on our 100% guarantee, please visit [www.novusbio.com/guarantee](http://www.novusbio.com/guarantee)

Earn gift cards/discounts by submitting a review: [www.novusbio.com/reviews/submit/362-0010](http://www.novusbio.com/reviews/submit/362-0010)

Earn gift cards/discounts by submitting a publication using this product:  
[www.novusbio.com/publications](http://www.novusbio.com/publications)

