Product Datasheet

NGFR/TNFRSF16/p75NTR Antibody 35180002-0.1mg

Unit Size: 0.1 mg

Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.

www.novusbio.com



technical@novusbio.com

Publications: 2

Protocols, Publications, Related Products, Reviews, Research Tools and Images at: www.novusbio.com/35180002

Updated 4/3/2021 v.20.1

Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications Submit a review at www.novusbio.com/reviews/destination/35180002



35180002-0.1mg

NGFR/TNFRSF16/p75NTR Antibody


Product Information	
Unit Size	0.1 mg
Concentration	Please see the vial label for concentration. If unlisted please contact technical services.
Storage	Store at 4C short term. Aliquot and store at -20C long term. Avoid freeze-thaw cycles.
Clonality	Polyclonal
Preservative	No Preservative
Isotype	lgG
Purity	Immunogen affinity purified
Buffer	20mM Potassium Phosphate (pH 7.0) and 0.15M NaCl

Product Description	
Host	Rabbit
Gene ID	4804
Gene Symbol	NGFR
Species	Human
Reactivity Notes	Human.
Specificity/Sensitivity	This product is specific for Human p75NTR/NGF
Immunogen	This antibody is specific for the C Terminus Region of the target protein.
Notes	Manufactured by SDIX's proprietary Genomic Antibody Technology™. GAT FAQs.

Product Application Details	
Applications	Western Blot, ELISA
Recommended Dilutions	Western Blot 1:5000-1:20000, ELISA 1:100-1:2000
Application Notes	This antibody is useful in ELISA and Western Blot.

Images

ELISA: NGF R/TNFRSF16/p75NTR Antibody [35180002]

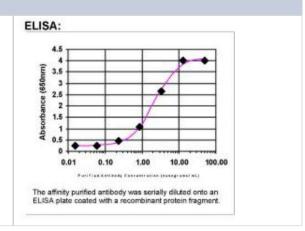

Publications

Garrison C, Lastwika K, Zhang Y et al. Proteomic Analysis, Immune Dysregulation, and Pathway Interconnections With Obesity J Proteome Res. 2017-01-06 [PMID: 27769113] (MiAr)

Details:

Analysis is performed on plasma proteomic data to identify how obesity can alter pathways and to highlight the risk factor for disease in subjects with a high body mass index.

Rho JH, Lampe PD. High-throughput screening for native autoantigen-autoantibody complexes using antibody microarrays J Proteome Res. 2013-05-03 [PMID: 23541305] (MiAr)

Details:

A novel method using antibody microarrays is used to detect autoantibody-antigen complexes that can potentially be useful for detection and characterization of diseases.





Novus Biologicals USA

10730 E. Briarwood Avenue Centennial, CO 80112 USA

Phone: 303.730.1950 Toll Free: 1.888.506.6887

Fax: 303.730.1966

nb-customerservice@bio-techne.com

Bio-Techne Ltd

19 Barton Lane Abingdon Science Park Abingdon, OX14 3NB, United Kingdom Phone: (44) (0) 1235 529449

Free Phone: 0800 37 34 15 Fax: (44) (0) 1235 533420 info.EMEA@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave Toronto, ON M8Z 4E6 Canada

Phone: 905.827.6400 Toll Free: 855.668.8722 Fax: 905.827.6402

canada.inquires@bio-techne.com

General Contact Information

www.novusbio.com Technical Support: nb-technical@bio-

techne.com
Orders: nb-customerservice@bio-techne.com

General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/35180002

Earn gift cards/discounts by submitting a publication using this product: www.novusbio.com/publications

