

# Product Datasheet

## Lightning-Link (R) Rapid Alexa Fluor 647 Antibody Labeling Kit 336-0030

Unit Size: 3 x 10ug Reaction

Store at -20C.

[www.novusbio.com](http://www.novusbio.com)



[technical@novusbio.com](mailto:technical@novusbio.com)

### Publications: 1

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:  
[www.novusbio.com/336-0030](http://www.novusbio.com/336-0030)

Updated 10/23/2024 v.20.1

Earn rewards for product  
reviews and publications.

Submit a publication at [www.novusbio.com/publications](http://www.novusbio.com/publications)

Submit a review at [www.novusbio.com/reviews/destination/336-0030](http://www.novusbio.com/reviews/destination/336-0030)




**336-0030****Lightning-Link (R) Rapid Alexa Fluor 647 Antibody Labeling Kit**

<b>Product Information</b>	
<b>Unit Size</b>	3 x 10ug Reaction
<b>Concentration</b>	Concentration is not relevant for this product. Please see the protocols for proper use of this product.
<b>Storage</b>	Store at -20C.
<b>Conjugate</b>	Alexa Fluor 647
<b>Product Description</b>	
<b>Description</b>	<p>Lightning-Link Rapid Antibody Labeling kits enable direct labeling of antibodies, proteins, peptides or other biomolecules (see protocol for further information) and allow conjugations to be set up in seconds, and used within minutes (less than 20 minutes).</p> <p>Our Lightning-Link Rapid technology can be used to label small quantities of antibody with 100% recovery. The researcher simply pipettes the antibody or other biomolecule into a vial of lyophilized mixture containing the label of interest and incubates the reaction for <b>only 15 minutes</b> (for more details please watch the video below). As there are no desalting or dialysis steps that commonly interrupt traditional protein labeling procedures.</p> <p><b>Features Benefits</b> Quick and easy to use Save time, no special knowledge required No separation steps 100% recovery - no antibody/protein loss Can be used in a wide range of applications Flexible Freeze dried Ships at ambient temperature, long shelf-life Fully scalable (10 ug to 1 g or more) Easy transfer from R&amp;D to manufacturing Stringently QC tested Consistent high quality, excellent batch-to-batch reproducibility Large number of labels available Experimental flexibility Reliable: nearly 300 references Successfully used in many fields of research</p> <p>Learn more about Lightning-Link™ Conjugation Kits by reading <a href="#">FAQs</a></p> <p>For more information please check out these useful links!  <a href="#">Antibody Labeling Guide</a>  <a href="#">Antibody Conjugation Illustrated Assay</a></p>
<b>Kit Components</b>	3 glass vials of Lightning-Link® Rapid mix, 1 vial of LL Rapid Modifier reagent, 1 vial of LL Rapid Quencher reagent
<b>Notes</b>	<p>This product is manufactured by Abcam and distributed by Novus Biologicals.</p> <p>This product is for research use only and is not approved for use in humans or in clinical diagnosis. This product is guaranteed for 1 year from date of receipt and this statement overrides any mentioned guarantee period on the limitations section of this products datasheet. Please contact <a href="mailto:technical@novusbio.com">technical@novusbio.com</a> with questions.</p>
<b>Product Application Details</b>	
<b>Application Notes</b>	By circumventing the desalting or dialysis steps that commonly interrupt traditional antibody conjugation procedures, LightningLink technology can be used to label both small (e.g. 10 ug) and large quantities of primary antibodies with ease. Batch-to-batch variation upon scale up is minimal as the process is so simple, and recoveries are always 100%.



## Images

Lightning-Link (R) Rapid Alexa Fluor 647 Antibody Labeling Kit [336-0030] - Scans of the Absorbance and Fluorescence Emission of a Lightning-Link® Rapid Alexa Fluor® 647 Conjugate

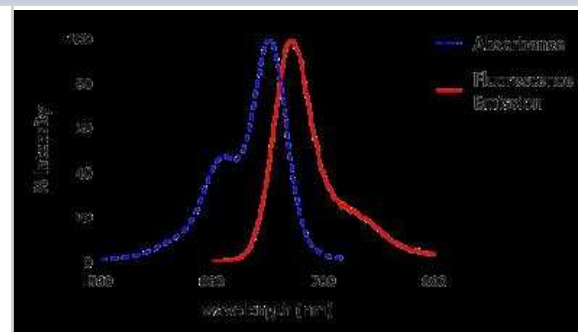

## Publications

Tanaka T, Warner BM, Michael DG Et al. LAMP3 inhibits autophagy and contributes to cell death by lysosomal membrane permeabilization Autophagy 2021-11-22 [PMID: 34802379]



### **Novus Biologicals USA**

10730 E. Briarwood Avenue  
Centennial, CO 80112  
USA  
Phone: 303.730.1950  
Toll Free: 1.888.506.6887  
Fax: 303.730.1966  
nb-customerservice@bio-techne.com

### **Bio-Techne Canada**

21 Canmotor Ave  
Toronto, ON M8Z 4E6  
Canada  
Phone: 905.827.6400  
Toll Free: 855.668.8722  
Fax: 905.827.6402  
canada.inquires@bio-techne.com

### **Bio-Techne Ltd**

19 Barton Lane  
Abingdon Science Park  
Abingdon, OX14 3NB, United Kingdom  
Phone: (44) (0) 1235 529449  
Free Phone: 0800 37 34 15  
Fax: (44) (0) 1235 533420  
info.EMEA@bio-techne.com

### **General Contact Information**

[www.novusbio.com](http://www.novusbio.com)  
Technical Support: nb-technical@bio-techne.com  
Orders: nb-customerservice@bio-techne.com  
General: novus@novusbio.com

### **Limitations**

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Kits are guaranteed for 6 months from date of receipt.

For more information on our 100% guarantee, please visit [www.novusbio.com/guarantee](http://www.novusbio.com/guarantee)

Earn gift cards/discounts by submitting a review: [www.novusbio.com/reviews/submit/336-0030](http://www.novusbio.com/reviews/submit/336-0030)

Earn gift cards/discounts by submitting a publication using this product:  
[www.novusbio.com/publications](http://www.novusbio.com/publications)

