

Product Datasheet

ABO, Blood Group A Antigen Antibody (HE-14) [Alexa Fluor® 350] NBP2-47868AF350

Unit Size: 0.1 ml

Store at 4C in the dark.

www.novusbio.com



technical@novusbio.com


Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/NBP2-47868AF350

Updated 9/10/2023 v.20.1

**Earn rewards for product
reviews and publications.**

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/NBP2-47868AF350

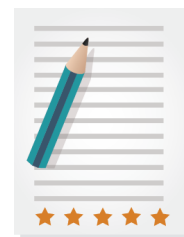

NBP2-47868AF350

ABO, Blood Group A Antigen Antibody (HE-14) [Alexa Fluor® 350]

Product Information	
Unit Size	0.1 ml
Concentration	Please see the vial label for concentration. If unlisted please contact technical services.
Storage	Store at 4C in the dark.
Clonality	Monoclonal
Clone	HE-14
Preservative	0.05% Sodium Azide
Isotype	IgM Kappa
Conjugate	Alexa Fluor 350
Purity	Protein A or G purified
Buffer	50mM Sodium Borate
Product Description	
Host	Mouse
Gene ID	28
Gene Symbol	ABO
Species	Human
Specificity/Sensitivity	This antibody recognizes human blood group A (monofucosyl and difucosyl A antigens with chain types 1, 2, 3, 4, 5, 6) and Forssmann antigen. It is also reactive with the immuno-dominant A trisaccharide. Blood group antigen expression in human colon cancer was studied by means of two monoclonal antibodies of broad anti-A (HE-14) and anti-type 3 and type 4 chain-based A and H (HE-10) specificity. These antigens were proved to re-appear in tumors of the distal colon, the HE-10 antibody reacting more frequently (9 out of 12 samples) than HE-14 (5 out of 12 samples) and frequently with supra-nuclear staining of the cytoplasm probably in those places of the Golgi apparatus where carbohydrate antigens are synthesized. This staining pattern is characteristic of HE-10 in normal colonic mucosa as well. With HE-14, staining was often absent in less differentiated tumors, while HE-10 did react in such tumors. In some cases, these two antibodies gave different staining patterns in parallel sections from the same tissue sample, primarily at the cellular level. Three out of 12 cases showed blood group antigen expression in the mucosa of the distal colon adjacent to the tumor only when HE-10 Ab was used.
Immunogen	Mixture of erythrocytes of blood group A1 and glyco protein fraction isolated from the saliva of secretors with blood group A

Notes	Alexa Fluor (R) products are provided under an intellectual property license from Life Technologies Corporation. The purchase of this product conveys to the buyer the non-transferable right to use the purchased product and components of the product only in research conducted by the buyer (whether the buyer is an academic or for-profit entity). The sale of this product is expressly conditioned on the buyer not using the product or its components, or any materials made using the product or its components, in any activity to generate revenue, which may include, but is not limited to use of the product or its components: (i) in manufacturing; (ii) to provide a service, information, or data in return for payment; (iii) for therapeutic, diagnostic or prophylactic purposes; or (iv) for resale, regardless of whether they are resold for use in research. For information on purchasing a license to this product for purposes other than as described above, contact Life Technologies Corporation, 5791 Van Allen Way, Carlsbad, CA 92008 USA or outlicensing@lifetech.com. This conjugate is made on demand. Actual recovery may vary from the stated volume of this product. The volume will be greater than or equal to the unit size stated on the datasheet.
--------------	--

Product Application Details	
Applications	Immunohistochemistry, Immunohistochemistry-Paraffin, Immunofluorescence
Recommended Dilutions	Immunohistochemistry, Immunohistochemistry-Paraffin, Immunofluorescence
Application Notes	Optimal dilution of this antibody should be experimentally determined.



Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112
USA

Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
novus@novusbio.com

Bio-Techne Canada

21 Canmotor Ave
Toronto, ON M8Z 4E6
Canada

Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB, United Kingdom
Phone: (44) (0) 1235 529449
Free Phone: 0800 37 34 15
Fax: (44) (0) 1235 533420
info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com
Technical Support: technical@novusbio.com
Orders: orders@novusbio.com
General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis.
Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP2-47868AF350

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

