

Product Datasheet

Cytokeratin 19 Antibody (KRT19/799 + KRT19/800)

[PE/Cy5.5]

NBP2-47948PECY55

Unit Size: 0.1 ml

Store at 4C in the dark. Do not freeze.

www.novusbio.com



technical@novusbio.com


Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/NBP2-47948PECY55

Updated 9/10/2023 v.20.1

Earn rewards for product
reviews and publications.

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/NBP2-47948PECY55

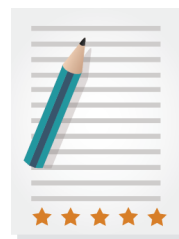

NBP2-47948PECY55

Cytokeratin 19 Antibody (KRT19/799 + KRT19/800) [PE/Cy5.5]

Product Information	
Unit Size	0.1 ml
Concentration	Please see the vial label for concentration. If unlisted please contact technical services.
Storage	Store at 4C in the dark. Do not freeze.
Clonality	Monoclonal
Clone	KRT19/799 + KRT19/800
Preservative	0.05% Sodium Azide
Isotype	IgG2a Lambda/IgG1 Kappa
Conjugate	PE/Cy5.5
Purity	Protein A or G purified
Buffer	PBS
Product Description	
Host	Mouse
Gene ID	3880
Gene Symbol	KRT19
Species	Human, Mouse
Marker	Pancreatic Stem Cell Marker
Specificity/Sensitivity	Recognizes a protein of 40kDa, identified as cytokeratin-19 (CK19), which is expressed in sweat gland, mammary gland ductal and secretory cells, bile ducts, gastrointestinal tract, bladder urothelium, oral epithelia, esophagus, and ectocervical epithelium. Anti-CK19 reacts with a wide variety of epithelial malignancies including adenocarcinomas of the colon, stomach, pancreas, biliary tract, liver, and breast. Perhaps the most useful application is the identification of thyroid carcinoma of the papillary type, although 50%-60% of follicular carcinomas are also labeled. Anti-CK19 is a useful marker for detection of tumor cells in lymph nodes, peripheral blood, bone marrow and breast cancer.
Immunogen	Recombinant full-length human Cytokeratin 19 protein (Uniprot: P08727)
Product Application Details	
Applications	Flow Cytometry
Recommended Dilutions	Flow Cytometry
Application Notes	Optimal dilution of this antibody should be experimentally determined. For optimal results using our Tandem dyes, please avoid prolonged exposure to light or extreme temperature fluctuations. These can lead to irreversible degradation or decoupling. When staining intracellular targets, specific attention to the fixation and permeabilization steps in your flow protocol may be required. Please contact our technical support team at technical@novusbio.com if you have any questions.





Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112
USA

Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
novus@novusbio.com

Bio-Techne Canada

21 Canmotor Ave
Toronto, ON M8Z 4E6
Canada

Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB, United Kingdom
Phone: (44) (0) 1235 529449
Free Phone: 0800 37 34 15
Fax: (44) (0) 1235 533420
info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com
Technical Support: technical@novusbio.com
Orders: orders@novusbio.com
General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP2-47948PECY55

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

