# **Product Datasheet**

# Bcl-6 Antibody (BCL6/1982) - Azide and BSA Free NBP2-75749

Unit Size: 100 ug

Store at -20 to -80C. Avoid freeze-thaw cycles.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at: www.novusbio.com/NBP2-75749

Updated 7/16/2024 v.20.1

Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications
Submit a review at www.novusbio.com/reviews/destination/NBP2-75749



# NBP2-75749


Bcl-6 Antibody (BCL6/1982) - Azide and BSA Free

Bcl-6 Antibody (BCL6/1982) - Azide and BSA Free	
Product Information	
Unit Size	100 ug
Concentration	1 mg/ml
Storage	Store at -20 to -80C. Avoid freeze-thaw cycles.
Clonality	Monoclonal
Clone	BCL6/1982
Preservative	No Preservative
Isotype	IgG1 Kappa
Purity	Protein A or G purified
Buffer	10 mM PBS
Target Molecular Weight	95 kDa
Product Description	
Description	1.0 mg/ml of antibody purified from Bioreactor Concentrate by Protein A/G. Prepared in 10mM PBS WITHOUT BSA & azide. Also available at 200 ug/ml WITH BSA & azide (NBP2-75748).  Antibody with azide - store at 2 to 8C. Antibody without azide - store at -20 to -80C.
Host	Mouse
Gene ID	604
Gene Symbol	BCL6
Species	Human
Marker	Follicular Lymphoma Marker
Specificity/Sensitivity	Recognizes a protein of 95kDa, which is identified as Bcl-6. Antibody to bcl-6 is helpful in a number of diagnostic settings: (1) In the differential diagnosis of small B-cell lymphoma. Follicular lymphoma will show bcl-6 (and CD10) positivity whereas other small B-cell lymphomas are usually negative. (2) Bcl-6 is an important prognostic marker in diffuse large B-cell lymphomas (DLBCL), where CD10, bcl-6 and MUM1/IRF4 are used to identify germinal center and activated B-cell phenotypes. (3) Bcl-6 can be valuable in distinguishing classical Hodgkin lymphoma from nodular lymphocyte predominant Hodgkin lymphoma (NLPHL). The Reed-Sternberg cells of classical Hodgkin lymphoma are bcl-6 negative whereas the large ('LH') cells of NLPHL are bcl-6 positive. In contrast, anti-Bcl-6 rarely stains mantle-cell lymphoma and MALT lymphoma.
Immunogen	Recombinant human bcl-6 protein fragment (around aa256-389) (exact sequence is proprietary) (Uniprot: P41182)
Product Application Details	
Applications	ELISA, Immunohistochemistry, Immunohistochemistry-Paraffin, Protein Array
Recommended Dilutions	ELISA 2 - 4 ug/mL, Immunohistochemistry 1 - 2 ug/mL, Immunohistochemistry-Paraffin 1 - 2 ug/mL, Protein Array
Application Notes	ELISA: Use Ab at 2-4ug/ml for coating. Order Ab without BSA. Immunohistochemistry (Formalin-fixed): 1-2ug/ml for 30 minutes at RT. Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95C followed by cooling at RT for 20 minutes. Optimal dilution for a specific application should be determined.




#### **Images**

Immunohistochemistry-Paraffin: Bcl-6 Antibody (BCL6/1982) - Azide and BSA Free [NBP2-75749] - Formalin-fixed, paraffin-embedded human Tonsil stained with BCL-6 Mouse Monoclonal Antibody (BCL6/1982).

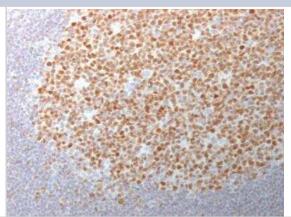

Protein Array: Bcl-6 Antibody (BCL6/1982) - Azide and BSA Free [NBP2-75749] - Z- and S- Score: The Z-score represents the strength of a signal that a monoclonal antibody (mAb) (in combination with a fluorescently-tagged anti-IgG secondary antibody) produces when binding to a particular protein on the HuProtTM array. Z-scores are described in units of standard deviations (SD's) above the mean value of all signals generated on that array. If targets on HuProtTM are arranged in descending order of the Z-score, the S-score is the difference (also in units of SD's) between the Z-score. S-score therefore represents the relative target specificity of a Monoclonal Antibody to its intended target. A mAb is considered to specific to its intended target, if the mAb has an S-score of at least 2.5. For example, if a mAb binds to protein X with a Z-score of 43 and to protein Y with a Z-score of 14, then the S-score for the binding of that mAb to protein X is equal to 29.

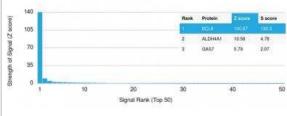



### **Novus Biologicals USA**

10730 E. Briarwood Avenue Centennial, CO 80112

USA

Phone: 303.730.1950 Toll Free: 1.888.506.6887

Fax: 303.730.1966

nb-customerservice@bio-techne.com

#### **Bio-Techne Canada**

21 Canmotor Ave Toronto, ON M8Z 4E6

Canada

Phone: 905.827.6400 Toll Free: 855.668.8722 Fax: 905.827.6402

Fax: 905.827.6402

canada.inquires@bio-techne.com

#### **Bio-Techne Ltd**

19 Barton Lane Abingdon Science Park Abingdon, OX14 3NB, United Kingdom Phone: (44) (0) 1235 529449

Free Phone: 0800 37 34 15 Fax: (44) (0) 1235 533420 info.EMEA@bio-techne.com

#### **General Contact Information**

www.novusbio.com

Technical Support: nb-technical@bio-

techne.com

Orders: nb-customerservice@bio-techne.com

General: novus@novusbio.com

# **Products Related to NBP2-75749**

HAF007 Goat anti-Mouse IgG Secondary Antibody [HRP]

NB720-B Rabbit anti-Mouse IgG (H+L) Secondary Antibody [Biotin]

NBP1-43319-0.5mg Mouse IgG1 Kappa Isotype Control (P3.6.2.8.1)
H00000604-Q01-10ug Recombinant Human Bcl-6 GST (N-Term) Protein

#### Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP2-75749

Earn gift cards/discounts by submitting a publication using this product: www.novusbio.com/publications

