

Product Datasheet

GITR Ligand/TNFSF18 Antibody (109101) [PE] FAB6941P-0.5ml

Unit Size: 0.5 ml

Store at 4C in the dark.

www.novusbio.com



technical@novusbio.com

Publications: 2

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/FAB6941P

Updated 3/10/2025 v.20.1

Earn rewards for product
reviews and publications.

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/FAB6941P



FAB6941P-0.5ml

GITR Ligand/TNFSF18 Antibody (109101) [PE]


Product Information	
Unit Size	0.5 ml
Concentration	Please see the vial label for concentration. If unlisted please contact technical services.
Storage	Store at 4C in the dark.
Clonality	Monoclonal
Clone	109101
Preservative	0.05% Sodium Azide
Isotype	IgG1
Conjugate	PE
Purity	Protein A or G purified from hybridoma culture supernatant
Buffer	PBS

Product Description	
Host	Mouse
Gene ID	8995
Gene Symbol	TNFSF18
Species	Human
Specificity/Sensitivity	Detects human GITR Ligand/TNFSF18 in flow cytometry.
Immunogen	<i>S. frugiperda</i> insect ovarian cell line Sf 21-derived recombinant human GITR Ligand/TNFSF18 Glu52-Ser177 Accession # Q9UNG2

Product Application Details	
Applications	Flow Cytometry
Recommended Dilutions	Flow Cytometry 10 uL/10 ⁶ cells
Application Notes	Optimal dilution of this antibody should be experimentally determined.

Images

HUVEC human umbilical vein endothelial cells were stained with Mouse Anti-Human GITR Ligand/TNFSF18 PE-conjugated Monoclonal Antibody (Catalog # FAB6941P, filled histogram) or isotype control antibody (Catalog # IC002P, open histogram). View our protocol for [Staining Membrane-associated Proteins](#).

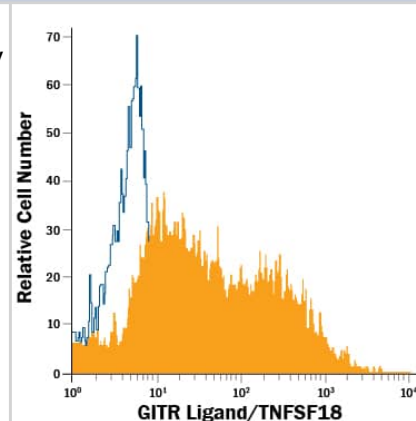

Publications

Wang C, Yi T, Qin L, Maldonado R, von Andrian U, Kulkarni S, Tellides G, Pober J Rapamycin-treated human endothelial cells preferentially activate allogeneic regulatory T cells. J Clin Invest, 2013-03-08;123(4):1677-93. 2013-03-08 [PMID: 23478407] (Flow Cytometry, Human)

M Toda, M Ueno, J Yamada, A Hiraga, K Asada, J Hamuro, C Sotozono, S Kinoshita Quiescent innate and adaptive immune responses maintain the long-term integrity of corneal endothelium reconstituted through allogeneic cell injection therapy Scientific Reports, 2022-10-27;12(1):18072. 2022-10-27 [PMID: 36302875] (Flow Cytometry, Human)



Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112

USA

Phone: 303.730.1950

Toll Free: 1.888.506.6887

Fax: 303.730.1966

nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave
Toronto, ON M8Z 4E6

Canada

Phone: 905.827.6400

Toll Free: 855.668.8722

Fax: 905.827.6402

canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane

Abingdon Science Park

Abingdon, OX14 3NB, United Kingdom

Phone: (44) (0) 1235 529449

Free Phone: 0800 37 34 15

Fax: (44) (0) 1235 533420

info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com

Technical Support: nb-technical@bio-techne.com

Orders: nb-customerservice@bio-techne.com

General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/FAB6941P

Earn gift cards/discounts by submitting a publication using this product:

www.novusbio.com/publications



