Product Datasheet

Heart Pericardium Tissue Slides (Lupus)- Frozen NBP2-77660

Unit Size: 5 Slides

Store at -70C. Avoid freeze-thaw cycles.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at: www.novusbio.com/NBP2-77660

Updated 1/9/2024 v.20.1

Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications Submit a review at www.novusbio.com/reviews/destination/NBP2-77660



NBP2-77660

Heart Pericardium Tissue Slides (Lupus)- Frozen


| Product Information | |
|---------------------|--|
| Unit Size | 5 Slides |
| Concentration | Concentration is not relevant for this product. Please see the protocols for proper use of this product. |
| Storage | Store at -70C. Avoid freeze-thaw cycles. |
| Draduct Description | |

| Storage | Store at -70C. Avoid freeze-thaw cycles. | |
|-----------------------------|--|--|
| Product Description | | |
| Description | Frozen tissue slides are ideal for rapid cellular localization of RNA and protein expression, especially for antibodies that don't recognize paraffin embedded tissue. Tissues were snap frozen in liquid nitrogen immediately after excision, and embedded in OCT. The tissue sections are 5-10um in thickness, mounted on positively charged glass slides, and fixed by cold acetone. At least one of the tissue slides from each lot was stained with H & E for quality assurance. | |
| Species | Human | |
| Notes | Donor information available upon request | |
| Lysate Type | Tissue | |
| Lysate Tissue | Heart | |
| Lysate Tissue Condition | Lupus | |
| Product Application Details | | |

| Product Application Details | |
|------------------------------------|--|
| Applications | Hematoxylin and Eosin Stain, Immunohistochemistry, Immunohistochemistry-Frozen, In-situ Hybridization |
| Recommended Dilutions | Immunohistochemistry, Immunohistochemistry-Frozen, In-situ Hybridization, Hematoxylin and Eosin Stain |
| Application Notes | Each individual tissue spot should be >90% core retention and total core lost from 1 slide should be no more than 5% |

Images

Hematoxylin & Eosin Stain: Heart Pericardium Tissue Slides (Lupus) [NBP2-77660] - Tissue: Human Heart (Pericardium), Pathology: Lupus

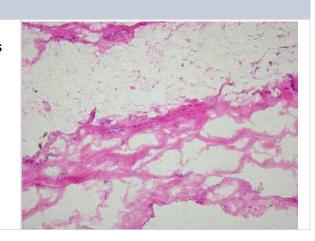



Novus Biologicals USA

10730 E. Briarwood Avenue Centennial, CO 80112 USA

Phone: 303.730.1950 Toll Free: 1.888.506.6887

Fax: 303.730.1966

nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave Toronto, ON M8Z 4E6 Canada

Phone: 905.827.6400 Toll Free: 855.668.8722 Fax: 905.827.6402

canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane Abingdon Science Park Abingdon, OX14 3NB, United Kingdom Phone: (44) (0) 1235 529449

Free Phone: 0800 37 34 15 Fax: (44) (0) 1235 533420 info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com Technical Support: nb-technical@bio-techne.com

Orders: nb-customerservice@bio-techne.com

General: novus@novusbio.com

Products Related to NBP2-77660

NB820-59226

Human Heart Pericardium Whole Tissue Lysate (Adult Whole Normal)

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Slides are guaranteed for 3 months from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP2-77660

Earn gift cards/discounts by submitting a publication using this product: www.novusbio.com/publications

