

Product Datasheet

DRAQ9 (TM) **NBP2-81128-200ul**

Unit Size: 200 ul

Store at -20C.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/NBP2-81128

Updated 2/24/2020 v.20.1

**Earn rewards for product
reviews and publications.**

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/NBP2-81128



NBP2-81128-200ul

DRAQ9 (TM)


Product Information	
Unit Size	200 ul
Concentration	Please see the protocols for proper use of this product. If no protocol is available, contact technical services for assistance.
Storage	Store at -20C.
Purity	>97%

Product Description	
Description	<p>Key Features of DRAQ9 (TM) : Far-red dye that stains membranous structures in the cytoplasm but not the nucleus Live-cell, permeant dye that can be used for live or fixed cells Track cells through mitotic cell division and monitor spheroid/organoid cultures in real-time Non-toxic and optimal for long-term imaging studies, up to several days Used for cell painting or mosaic analysis in high content phenotypic screening Spectrally compatible with GFP/FITC and UV-excited probes such as Hoechst</p>
Notes	Products are shipped at ambient temperature, but on receipt packs should be stored at -20C. Thaw and adjust temperature to above 20C before use.


Product Application Details	
Applications	Immunocytochemistry/Immunofluorescence, Live Imaging Microscopy, Fluorescence Imaging
Recommended Dilutions	Immunocytochemistry/Immunofluorescence 1:50 - 1:500, Live Imaging Microscopy 1:500, Fluorescence Imaging 1:50 - 1:500
Application Notes	DRAQ9 (TM) is supplied as a DMSO solution at a concentration of 1mM: the 200 ul size equates to 8x 24-multi-well plates and 5x 96-microtiter plates. The 1 ml size equates to 40x 24-multi-well plates, 25x 96-microtiter plates, and 3.3x 384-microtiter plates.

Images


Fluorescence Imaging: DRAQ9 (TM) [NBP2-81128] - Cells stained with .

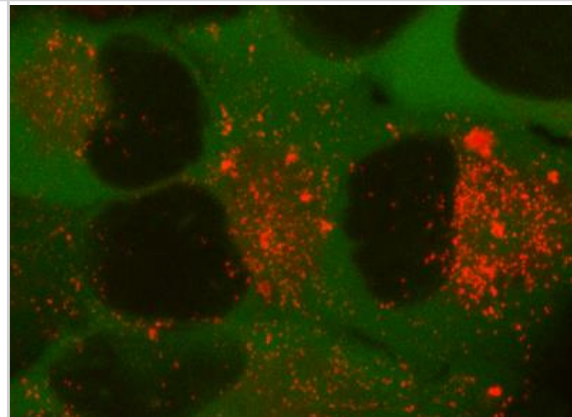


Live Imaging Microscopy: DRAQ9 (TM) [NBP2-81128] - DRAQ9 was diluted in EMEM culture media and applied to live HeLa cells at 2 μ M (1:500) for 30 minutes at room temperature and protected from light. Imaging was done immediately after staining without washing the cells.

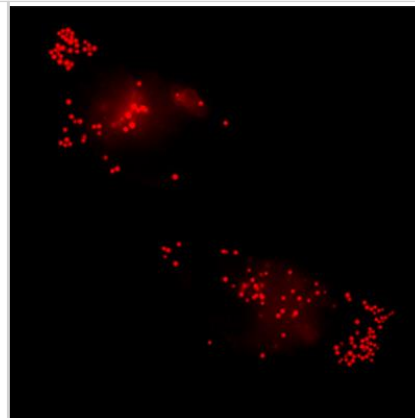


DRAQ9 (TM) [NBP2-81128] - Excitation and Emission Spectra for DRAQ9 (TM).

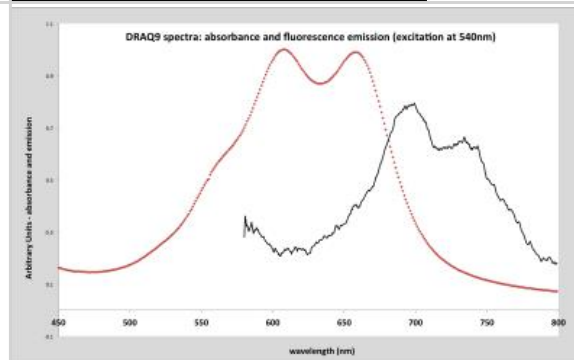

Live Imaging Microscopy: DRAQ9 (TM) [NBP2-81128] - labelling of live U2OS cells for 48 h at 2 μ M - showing a number of cells in mitosis.

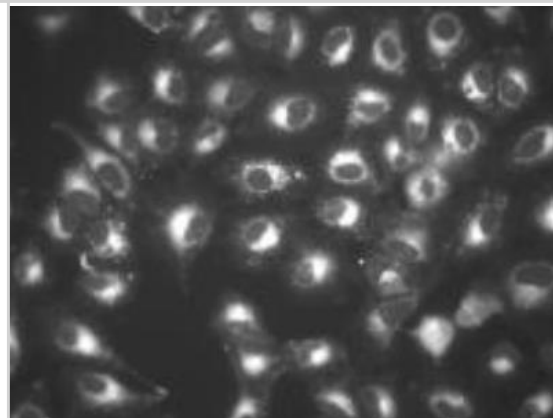

Fluorescence Imaging: DRAQ9 (TM) [NBP2-81128] - 6-day U2OS spheroid with continuous labeling with at 2 μ M.

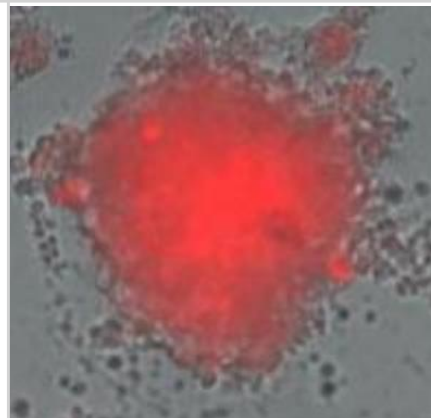



Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112
USA

Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
novus@novusbio.com

Bio-Techne Canada

21 Canmotor Ave
Toronto, ON M8Z 4E6
Canada

Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB, United Kingdom

Phone: (44) (0) 1235 529449
Free Phone: 0800 37 34 15
Fax: (44) (0) 1235 533420
info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com
Technical Support: technical@novusbio.com
Orders: orders@novusbio.com
General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Support products are guaranteed for 6 months from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP2-81128

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

