

Product Datasheet

Purple Hematoxylin Solution NBP3-12176

Unit Size: 100 ml

Store at room temperature.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/NBP3-12176

Updated 12/9/2021 v.20.1

Earn rewards for product
reviews and publications.

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/NBP3-12176



NBP3-12176


Purple Hematoxylin Solution

Product Information	
Unit Size	100 ml
Concentration	Please see the protocols for proper use of this product. If no protocol is available, contact technical services for assistance.
Storage	Store at room temperature.

Product Application Details	
Applications	Immunohistochemistry, Immunohistochemistry-Paraffin
Recommended Dilutions	Immunohistochemistry, Immunohistochemistry-Paraffin
Application Notes	<p>Preparation of Working Solutions</p> <p>Purple Hematoxylin is ready to use and does not require any preparation.</p> <p>Protocol Recommendations</p> <p>Counterstain</p> <ol style="list-style-type: none"> 1. Rinse slides with deionized or distilled water. 2. Place slides into Purple Hematoxylin for 30-60 seconds. 3. Rinse slides in deionized or distilled water until all free hematoxylin has been removed. 4. Bluing the slides in an alkaline solution is not required and is not recommended.

Images

Immunohistochemistry-Paraffin: Purple Hematoxylin Solution [NBP3-12176] - Formalin fixed paraffin embedded human colon stained with epithelial antigen antibody labeled with Purple Hematoxylin Solution (NBP3-12176)

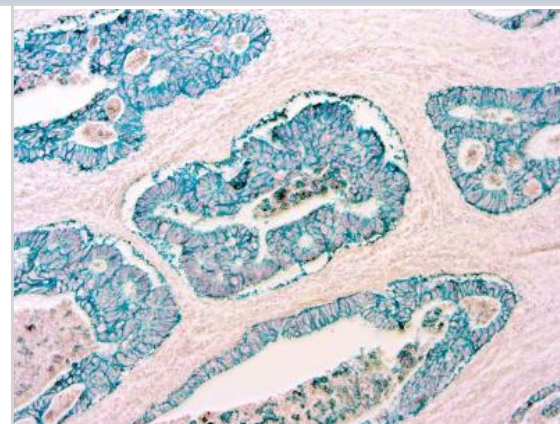



Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112
USA

Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
novus@novusbio.com

Bio-Techne Canada

21 Canmotor Ave
Toronto, ON M8Z 4E6
Canada

Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB, United Kingdom

Phone: (44) (0) 1235 529449
Free Phone: 0800 37 34 15
Fax: (44) (0) 1235 533420
info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com
Technical Support: technical@novusbio.com
Orders: orders@novusbio.com
General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Support products are guaranteed for 6 months from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP3-12176

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

