

Product Datasheet

Tissue Preservation Solution - HRP/AP assays NBP3-12178

Unit Size: 30 ml

Store at room temperature.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/NBP3-12178

Updated 12/9/2021 v.20.1

Earn rewards for product
reviews and publications.

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/NBP3-12178



NBP3-12178


Tissue Preservation Solution - HRP/AP assays

Product Information	
Unit Size	30 ml
Concentration	Please see the protocols for proper use of this product. If no protocol is available, contact technical services for assistance.
Storage	Store at room temperature.

Product Application Details	
Applications	Immunohistochemistry
Recommended Dilutions	Immunohistochemistry
Application Notes	<p>Preparation of Working Solutions Tissue Preservation Solution - HRP/AP assays is ready to use and does not require additional preparation.</p> <p>Protocol Recommendations</p> <ol style="list-style-type: none"> 1. After counterstaining, rinse slide with water. 2. Position the slide at an angle on a stack of paper towels to allow excess fluid to run off. 3. Lay slides flat and pipette or add a few drops of Tissue Preservation Solution - HRP/AP assays over the tissue, making sure you cover the tissue section. 4. Let tissue section soak for about 5 minutes. 5. Position the slide at an angle on a stack of paper towels to allow excess fluid to run off. 6. Let the fluid dry completely without further rinsing for about a minimum of 20 minutes or longer if needed. It is important that the slide is completely dry before moving on to the next step. 7. Mount and coverslip slide using any conventional mounting method.

Images

Tissue Preservation Solution - HRP/AP assays [NBP3-12178]

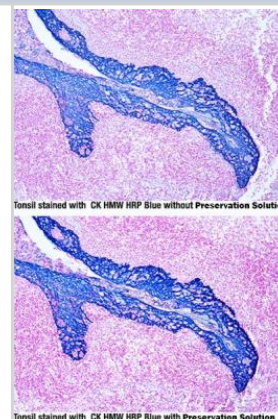



Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112
USA

Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
novus@novusbio.com

Bio-Techne Canada

21 Canmotor Ave
Toronto, ON M8Z 4E6
Canada

Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB, United Kingdom

Phone: (44) (0) 1235 529449
Free Phone: 0800 37 34 15
Fax: (44) (0) 1235 533420
info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com
Technical Support: technical@novusbio.com
Orders: orders@novusbio.com
General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Support products are guaranteed for 6 months from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP3-12178

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

