Product Datasheet

BAD-LAMP/LAMP5 Antibody (LAMP5/7642) - Azide and BSA Free NBP3-24264

Unit Size: 100 ug

Store at -20 to -80C. Avoid freeze-thaw cycles.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at: www.novusbio.com/NBP3-24264

Updated 7/16/2024 v.20.1

Earn rewards for product reviews and publications.

Submit a publication at www.novusbio.com/publications
Submit a review at www.novusbio.com/reviews/destination/NBP3-24264



NBP3-24264


BAD-LAMP/LAMP5 Antibody (LA	MP5/7642) - Azide and BSA Free
Product Information	
Unit Size	100 ug
Concentration	1 mg/ml
Storage	Store at -20 to -80C. Avoid freeze-thaw cycles.
Clonality	Monoclonal
Clone	LAMP5/7642
Preservative	No Preservative
Isotype	IgG2 Kappa
Purity	Protein A or G purified
Buffer	10mM PBS
Product Description	
Description	Antibody without azide - store at -20 to -80C.
Host	Mouse
Gene ID	24141
Gene Symbol	LAMP5
Species	Human
Marker	Late Endosomes Marker
Immunogen	Recombinant full-length human BAD-LAMP/LAMP5 protein
Product Application Details	
Applications	Western Blot, Immunohistochemistry-Paraffin
Recommended Dilutions	Western Blot 1-2ug/ml, Immunohistochemistry-Paraffin 1-2ug/ml
Application Notes	Immunohistochemistry-Paraffin: 1-2ug/ml for 30 minutes at RT. Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA pH 9.0 for 45 min at 95C followed by cooling at RT for 20 minutes

Applications	Western Blot, Immunohistochemistry-Paraffin
Recommended Dilutions	Western Blot 1-2ug/ml, Immunohistochemistry-Paraffin 1-2ug/ml
Application Notes	Immunohistochemistry-Paraffin: 1-2ug/ml for 30 minutes at RT. Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95C followed by cooling at RT for 20 minutes. Optimal dilution for a specific application should be determined.




Images

Formalin-fixed, paraffin-embedded human kidney stained with BAD-LAMP/LAMP5 antibody (LAMP5/7642) - Azide and BSA Free. Inset: PBS instead of primary antibody; secondary only negative control.

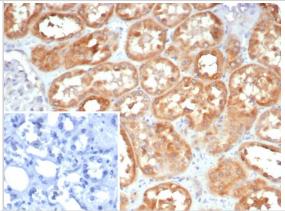

Western blot analysis of U937 cell lysate using BAD-LAMP/LAMP5 antibody (LAMP5/7642) - Azide and BSA Free.

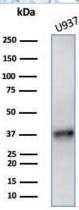



Novus Biologicals USA

10730 E. Briarwood Avenue Centennial, CO 80112

USA

Phone: 303.730.1950 Toll Free: 1.888.506.6887

Fax: 303.730.1966

nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave Toronto, ON M8Z 4E6

Canada

Phone: 905.827.6400 Toll Free: 855.668.8722 Fax: 905.827.6402

canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane Abingdon Science Park Abingdon, OX14 3NB, United Kingdom Phone: (44) (0) 1235 529449

Free Phone: 0800 37 34 15 Fax: (44) (0) 1235 533420 info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com

Technical Support: nb-technical@bio-

techne.com

Orders: nb-customerservice@bio-techne.com

General: novus@novusbio.com

Products Related to NBP3-24264

HAF007 Goat anti-Mouse IgG Secondary Antibody [HRP]

NB720-B Rabbit anti-Mouse IgG (H+L) Secondary Antibody [Biotin]

NBP1-84246PEP BAD-LAMP/LAMP5 Recombinant Protein Antigen

NBL1-08354 BAD-LAMP/LAMP5 Overexpression Lysate

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Primary Antibodies are guaranteed for 1 year from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP3-24264

Earn gift cards/discounts by submitting a publication using this product: www.novusbio.com/publications

