

Product Datasheet

Blu13 Prestained Protein Ladder NBP3-33174

Unit Size: 500 ul

Store at -20C. Avoid freeze-thaw cycles.

www.novusbio.com



technical@novusbio.com

Protocols, Publications, Related Products, Reviews, Research Tools and Images at:
www.novusbio.com/NBP3-33174

Updated 10/17/2024 v.20.1

Earn rewards for product
reviews and publications.

Submit a publication at www.novusbio.com/publications

Submit a review at www.novusbio.com/reviews/destination/NBP3-33174



NBP3-33174**Blu13 Prestained Protein Ladder****Product Information**

Unit Size	500 ul
Concentration	Please see the protocols for proper use of this product. If no protocol is available, contact technical services for assistance.
Storage	Store at -20C. Avoid freeze-thaw cycles.
Buffer	20 mM Trisphosphate pH 7.5 at 25C, 2% SDS, 0.2 mM Dithiothreitol, 3.6 M Urea, 15% (v/v) Glycerol.

Product Description


Description	The Blu13 Prestained Protein Ladder is a combination of 13 pre-stained proteins with molecular weights from 5 to 245 kDa. The 13 recombinant proteins are covalently coupled with blue chromophore, while 1 green band at 25 kDa and a red band at 75 kDa serve as reference bands. The Blu13 Prestained Protein Ladder keeps track of the size and separation of proteins during SDS-polyacrylamide gel electrophoresis, approximating protein size and validating Western transfer efficiency on PVDF, nylon, or nitrocellulose membranes. (13 pre-stained bands, 5-245 kDa)
--------------------	--

Product Application Details

Applications	Western Blot
Recommended Dilutions	Western Blot
Application Notes	3 ul or 5 ul per loading for clear visualization during electrophoresis on 15-well or 10-well mini-gel, respectively. 2.5 ul per well for general Western transferring.

Images

Western Blot: Blu13 Prestained Protein Ladder [NBP3-33174] - The Blu13 Prestained Protein Ladder is a combination of 13 pre-stained proteins with molecular weights from 5 to 245 kDa. The 13 recombinant proteins are covalently coupled with blue chromophore, while 1 green band at 25 kDa and a red band at 75 kDa serve as reference bands. The Blu13 Prestained Protein Ladder keeps track of the size and separation of proteins during SDS-polyacrylamide gel electrophoresis, approximating protein size and validating Western transfer efficiency on PVDF, nylon, or nitrocellulose membranes.

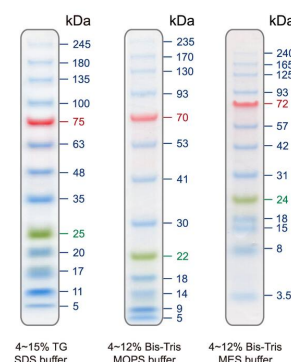



Novus Biologicals USA

10730 E. Briarwood Avenue
Centennial, CO 80112
USA
Phone: 303.730.1950
Toll Free: 1.888.506.6887
Fax: 303.730.1966
nb-customerservice@bio-techne.com

Bio-Techne Canada

21 Canmotor Ave
Toronto, ON M8Z 4E6
Canada
Phone: 905.827.6400
Toll Free: 855.668.8722
Fax: 905.827.6402
canada.inquires@bio-techne.com

Bio-Techne Ltd

19 Barton Lane
Abingdon Science Park
Abingdon, OX14 3NB, United Kingdom
Phone: (44) (0) 1235 529449
Free Phone: 0800 37 34 15
Fax: (44) (0) 1235 533420
info.EMEA@bio-techne.com

General Contact Information

www.novusbio.com
Technical Support: nb-technical@bio-
techne.com
Orders: nb-customerservice@bio-techne.com
General: novus@novusbio.com

Limitations

This product is for research use only and is not approved for use in humans or in clinical diagnosis. Support products are guaranteed for 6 months from date of receipt.

For more information on our 100% guarantee, please visit www.novusbio.com/guarantee

Earn gift cards/discounts by submitting a review: www.novusbio.com/reviews/submit/NBP3-33174

Earn gift cards/discounts by submitting a publication using this product:
www.novusbio.com/publications

